

MAMMOMAT 1000/3000/3000 Nova - Opdima®



Service

Service Instructions

© Siemens AG 2001

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

Register 5

Print No.: SPB7-230.061.08.03.02

Replaces: SPB7-230.061.08.02.02 and SPB7-230.061.06.02.02

English

Doc. Gen. Date: 03.01

65 69 334

Chapter	Page	Revision
All	All	03

1	Prerequisites	1 - 1
	General	1 - 1
	Training of customer support engineers.	1 - 1
	Documents required.	1 - 1
	Ultra 10 CD-ROMs required.	1 - 1
	Ultra 1 CD-ROMs required	1 - 1
	Meters and appliances required.	1 - 2
	Tools required	1 - 2
2	Functional description	2 - 1
	Ultra 10 Block diagram	2 - 1
	Ultra 1 Block diagram	2 - 2
	Biopsy controller.	2 - 3
	Biopsy unit control	2 - 3
	RS-232 switching	2 - 3
	CCD camera power supply	2 - 3
	LEDs.	2 - 3
	Workstation	2 - 4
	Application software	2 - 4
	CCD camera communication	2 - 5
	Control of MAMMOMAT and biopsy unit	2 - 5
	Monitor.	2 - 5
	Power supply workstation main unit.	2 - 5
	CCD camera.	2 - 5
	Biopsy unit.	2 - 5
3	Protective measures for CCD camera	3 - 1
4	Protective measures	4 - 1
5	Service mode	5 - 1
	General	5 - 1
	Advanced service	5 - 1
	Selection of mode	5 - 2
	Calibration of the biopsy unit	5 - 3
	Setting up needles.	5 - 6
	Editing needle values	5 - 7
	Backup of temporary storage media	5 - 8
	Use of Advanced service functions	5 - 9
	General	5 - 9
	Main menu.	5 - 9
	Country settings	5 - 11
	Test of units	5 - 12
	Software upgrade	5 - 13
	New software version	5 - 14
	Enabling the DICOM option	5 - 15

	Page
Instruction for obtaining the systems host ID and hostname5 - 15
Network settings5 - 17
Miscellaneous.5 - 18
Restoring data from MO disk5 - 20
CCD camera calibration and maintenance5 - 21
Grid table5 - 21
Non grid table.5 - 22
Disk cache settings5 - 23
Use of printer setup.5 - 24
Stand alone Opdima system5 - 24
Networked Opdima system5 - 25
Log administration5 - 27
Log inspection5 - 28
6 Removal and replacement of sub-assemblies	6 - 1
General	6 - 1
MO unit (no CD)	6 - 1
Removing the MO unit	6 - 1
Installing the MO unit	6 - 2
Tests and Adjustments	6 - 2
CD/MO unit	6 - 3
Removing the CD/MO unit	6 - 3
Installing the CD/MO	6 - 4
Tests and Adjustments	6 - 4
Biopsy controller	6 - 5
Removal of biopsy controller	6 - 5
Removal of biopsy controller cover	6 - 6
Replacement of components in the biopsy controller.	6 - 6
Ultra 10 Cables	6 - 8
Ultra 1 Cables	6 - 9
Biopsy unit	6 - 10
Ultra 10 Workstation	6 - 11
Connectors	6 - 11
Ultra 1 Workstation	6 - 12
Connectors	6 - 12
CCD camera	6 - 13
Exchange of workstation	6 - 14
Procedure.	6 - 14
Reinstallation of software.	6 - 16
Ultra 10 Pre installation	6 - 16
Ultra 1 Pre installation	6 - 20
Ultra 10/Ultra 1 Installation of Opdima ASW 2.1 (part No. 65 27 639)	6 - 22
Restore hostname and data	6 - 24
Final procedures	6 - 24
7 Messages	7 - 1
Workstation	7 - 1
Control panel	7 - 3

8	Fault isolation chart	8 - 1
	Fault isolation chart	8 - 1
	Explanations.	8 - 2
9	Troubleshooting guide	9 - 1
	General	9 - 1
	Biopsy unit.	9 - 2
	Biopsy unit not responding	9 - 2
	Problem with calibration	9 - 2
	Camera	9 - 3
	Camera not responding	9 - 3
	Cannot use camera	9 - 3
	Workstation	9 - 4
	If "Bogus file system" appears	9 - 4
	Problem with date in database mode	9 - 4
	No images displayed on monitor	9 - 4
	Tests.	9 - 4
	Ultra 1 If the workstation does not start after an unplanned power cut off.	9 - 4
	Image Quality	9 - 5
	Problems during calibration	9 - 5
	Quadrant difference on patient images	9 - 5
	Quadrant missing	9 - 5
	White line/dot	9 - 5
	MAMMOMAT	9 - 6
	Cannot use MAMMOMAT	9 - 6
	Network problem	9 - 7
	Printer	9 - 8
	Printer connected directly to Opdimia	9 - 8
	Printer connected to network	9 - 8
	MO unit / CD drive.	9 - 9
	All problems with MO disk	9 - 9
	If the message "Failed to store on MO disk" appears	9 - 9
	If the message "Cannot read disk Prepare disk?" appears	9 - 9
	Cannot communicate with MO unit	9 - 9
10	Measures after service	10 - 1
	Verifying the calibration of the biopsy unit.	10 - 1
	Check of Opdimia AEC	10 - 2
	Procedure	10 - 2
	Performance Criteria.	10 - 3
	Check of resolution	10 - 4
	Procedure	10 - 4
	Performance Criteria.	10 - 4
	Protective earth measurement	10 - 5
	Biopsy unit.	10 - 5

	Biopsy controller10 - 5
	CCD camera10 - 5
	Biopsy controller cable duct10 - 6
11	Changes to previous version	11 - 1
12	Appendix 1	A1 - 1
	Ultra 1 Installation of CD drive, part no 65 38 180	A1 - 1
	Ultra 1 Cables	A1 - 2
13	Appendix 2	A2 - 1
	Database log file	A2 - 1
14	Appendix 3	A3 - 1
	Customer specific data	A3 - 1
	Test protocol repetitiveness	A3 - 1
	Biopsy calculations	A3 - 1
	Test protocol CCD camera calibration	A3 - 2
	Test protocol image quality	A3 - 3
	AEC function	A3 - 3
	Resolution	A3 - 5
	Test protocol protective earth measurement	A3 - 6
	Biopsy unit	A3 - 6
	Biopsy controller	A3 - 6
	CCD camera	A3 - 6
	Biopsy controller cable duct	A3 - 6

General

Valid for Opdima® system ASW 2.1 on SUN workstation, part No. 64 30 453.

This document is valid for both Ultra 10 (SUN workstation with serial No. ≥ 1201) and Ultra 1 (SUN workstation with serial No. < 1201) when not stated differently.

Training of customer support engineers

Due to the technology used in this equipment, setup, service and maintenance may only be carried out by a customer support engineer who has attended a training workshop or has participated in at least one installation.

Documents required

- Supplement to the Instructions for Use MAMMOMAT 3000 - Opdima® (included in the Opdima® delivery)
- MAMMOMAT 1000/3000/3000 Nova - Opdima® Maintenance Instructions (included in the Opdima® delivery)
- MAMMOMAT 3000 - Opdima® Wiring Diagram (included in the Opdima® delivery)
- MAMMOMAT 1000/3000 Nova Wiring Diagram
- Technical Information TI 219
- **For Ultra 1 only:** MAMMOMAT 3000 Opdima® Service Instructions Workstation (Ultra 1) SPB7-230.061.02.02.02

Ultra 10 CD-ROMs required

- Ultra 10 Service Manual, part No. 65 27 670 (included in the Opdima® delivery)

Reading Ultra 10 Service Manual CD-ROM:

You may view the HTML document and watch the video clips on a PC with a CD drive running Windows 95 or later by using Netscape Navigator™, version 4.0 or later, or Internet Explorer, version 4.0 or later. Refer to the browser documentation about configuring your browser to view MPEG-1 files.

- Operating environment installation CD 02/00, part No. 65 27 704 or Operating environment installation CD 05/99, part No. 65 27 696 (for Solaris 2.5.1 HW 11/97, included in the Opdima® delivery)
- Solaris 2.5.1 HW 11/97 CD-ROM and CDE/Desktop 1.1 CD-ROM (both are included in a separate package called Media kit Solaris 2.5.1 HW 11/97, part No. 65 01 956, included in the workstation delivery)
- Opdima® ASW 2.1, part No. 65 27 639

Ultra 1 CD-ROMs required

- Opdima® ASW 1.3, part No. 64 32 665 (includes operating system)
- Opdima® ASW 2.1, part No. 65 27 639

Meters and appliances required


- Protective ground wire tester (44 15 899 RV090)
- Stereo calibration phantom (included in the Opdima® delivery) (part No. 64 30 701)
- Acrylic plastic (PMMA) plates, three plates measuring 150 mm x 150 mm x 20 mm and one plate measuring 150 mm x 150 mm x 10 mm, part No. 85 49 438 and 97 88 423 respectively
- Resolution phantom with at least 10 line pairs per mm
- 4,5 cm PMMA
- **For Ultra 10 only:** PC with a CD drive running Windows 95 or later by using Netscape Navigator™, version 4.0 or later, or Internet Explorer, version 4.0 or later (to be able to read Ultra 10 Service Manual)
- **If software upgrade:** External CD drive for Ultra 1 (part No. 64 30 453 serial No. 1109-1146)

Tools required

- Standard service tools
- PROM extractor for PLCC 32

Ultra 10 Block diagram

The block diagram below shows the function of the Opdima® system.

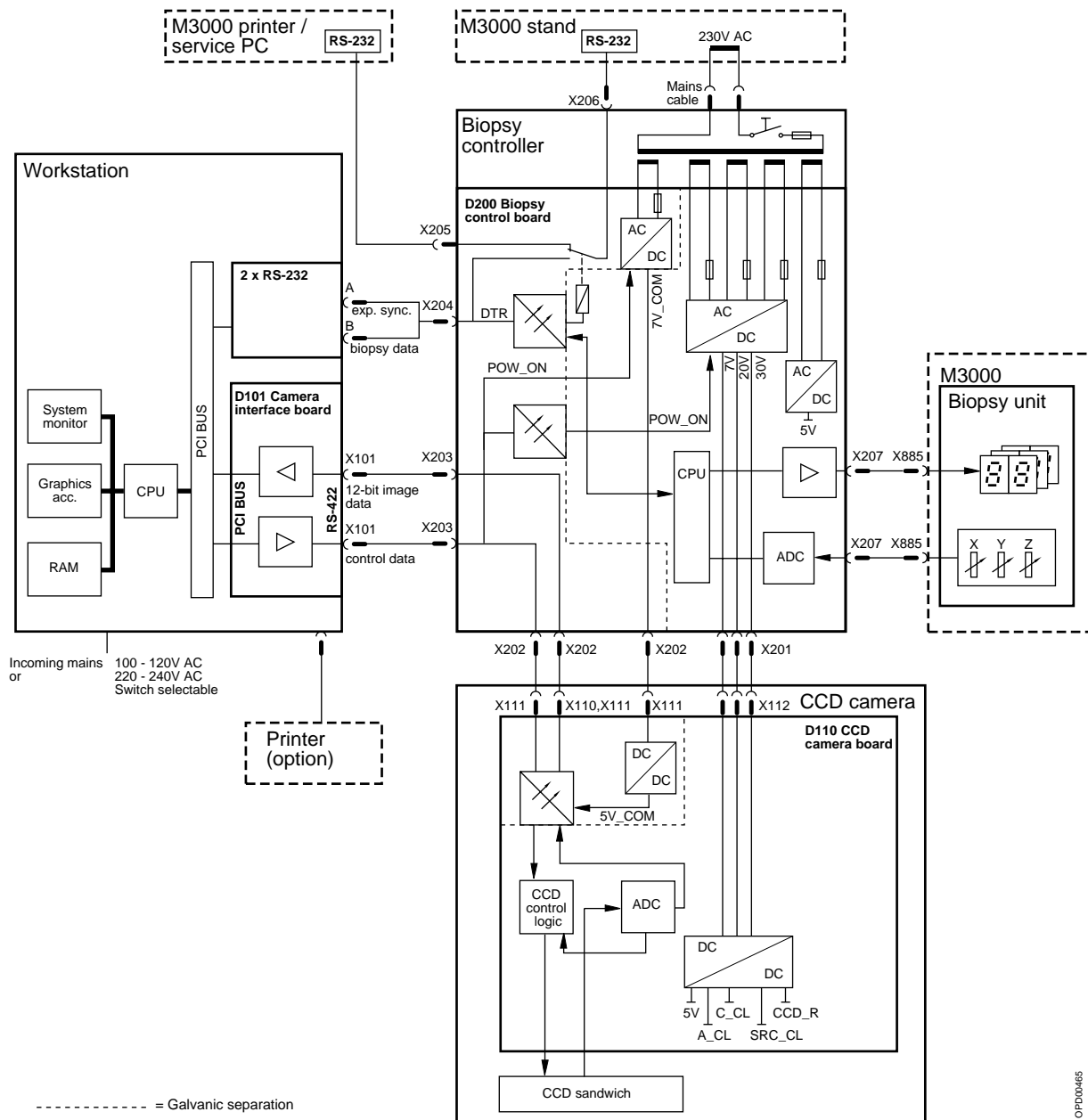

OPD00465

Fig. 1 Block diagram

Ultra 1 Block diagram

The block diagram below shows the function of the Opdima® system.

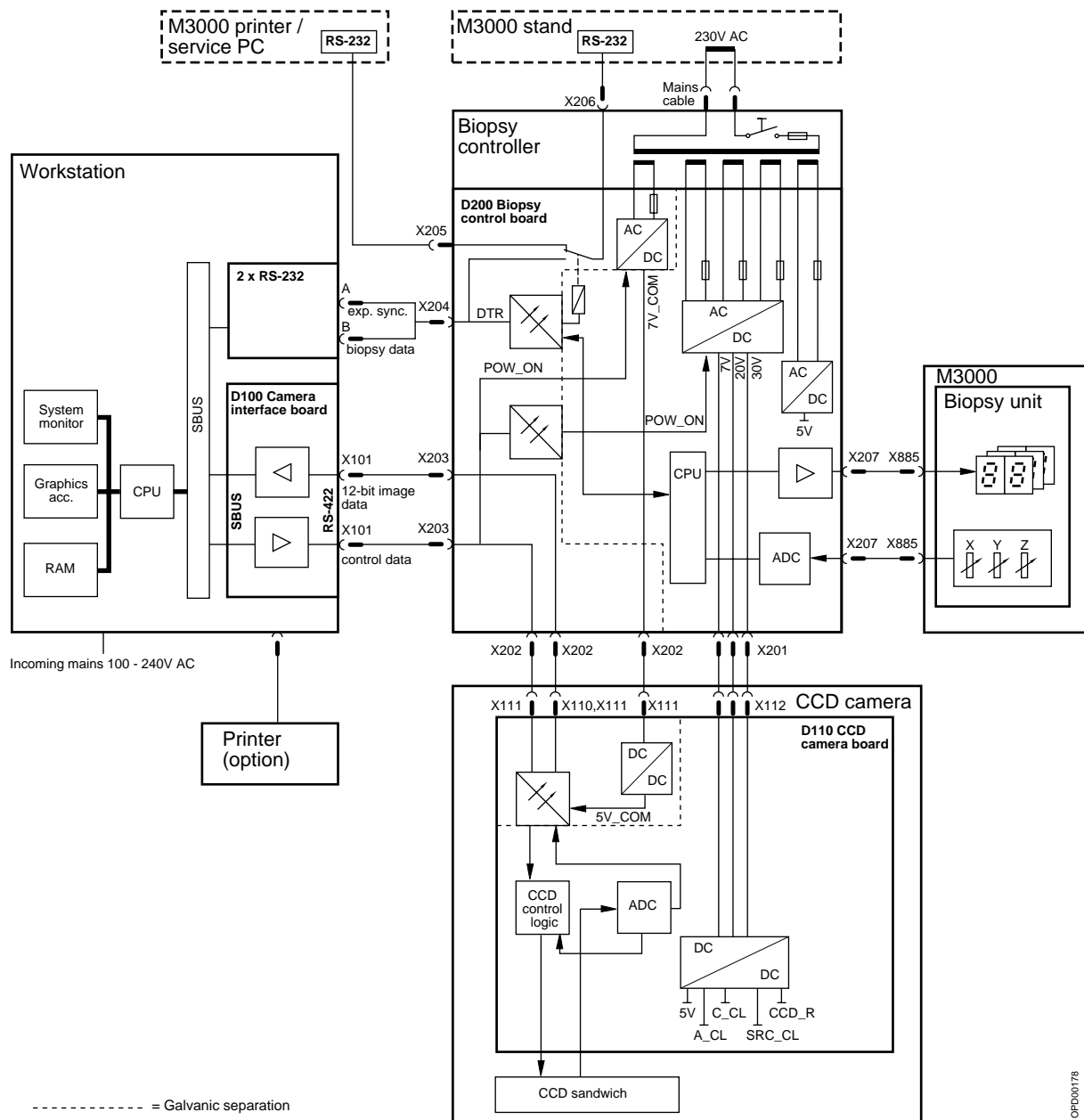


Fig. 2 Block diagram

OPD00178

Biopsy controller

The biopsy controller has three main functions:

- Biopsy unit control
- Switching of RS-232 from MAMMOMAT to either MAMMOMAT printer/service PC or workstation
- CCD camera power supply

Biopsy unit control

The biopsy functions are controlled by the micro controller. In the biopsy unit there are three potentiometers indicating the current needle position.

The values from the potentiometers are A/D converted and sent to the micro controller. To check the converted values from the A/D converter, the potentiometer values are converted in parallel by the micro controller.

RS-232 switching

To enable RS-232 communication between MAMMOMAT and either MAMMOMAT printer/service PC or workstation, a relay is used. The relay is controlled by the workstation.

CCD camera power supply

The biopsy controller converts the mains AC voltage to DC voltages. The DC voltages are used by the CCD camera.

LEDs

On the biopsy controller printed circuit board, D200, there are a number of LEDs with the following meaning:

- Indication of RAM error
- Indication of PROM error
- Indication of analog/digital converter
- Indication of NVM error
- Indication of TxD error
- Selection of RS-232 communication to MAMMOMAT printer/service PC or workstation (SELECT_WS)
- Supply voltages (e.g. 5V_COM, 30V_CC)

For more information on LEDs, see MAMMOMAT 1000/3000/3000 Nova - Opdima® Wiring Diagram.

Workstation

The workstation includes the following parts:

- Main unit, including temporary storage media i.e. hard disk
- Monitor
- Keyboard and mouse
- CD drive
 - For Ultra 1 only:** External CD drive
 - For Ultra 10 only:** Internal CD drive
- Solaris operating system
 - For Ultra 1 only:** Included on the ASW 1.3 CD-ROM
 - For Ultra 10 only:** Included in Media kit Solaris 2.5.1 HW 11/97
- External MO unit i.e. magneto-optical drive, used for the permanent storage media and backup
- Printer (option)

The workstation is used for:

- Running application software
- Communication with the CCD camera
- Control of MAMMOMAT and biopsy unit

Application software

The application software run by the workstation is the main communication interface with the user and units of the Opdima® system.

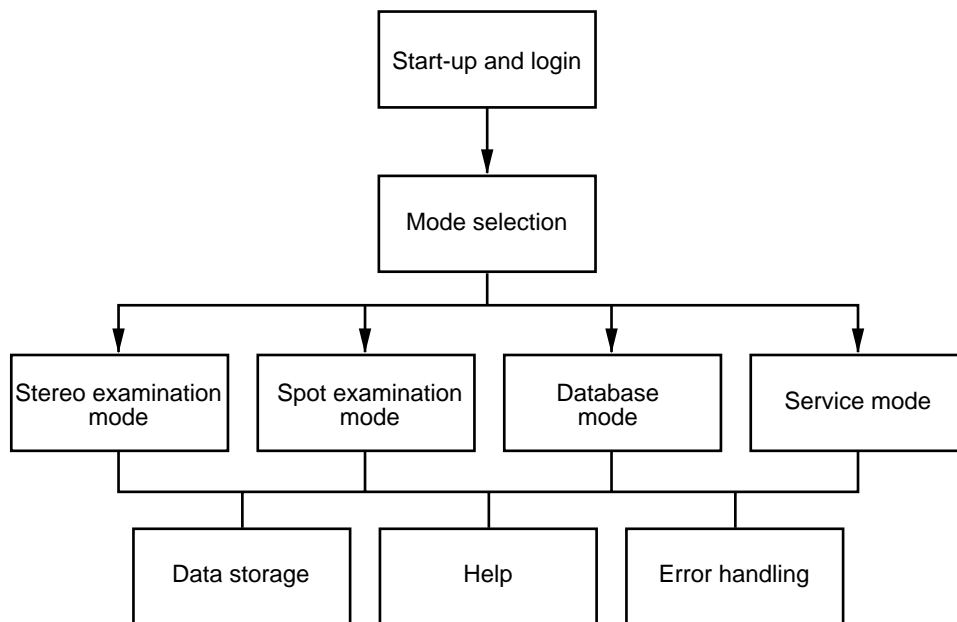


Fig. 3 Application software logic flow

CCD camera communication

The camera interface in the workstation main unit sends commands to the CCD camera from the workstation and receives the image from the camera. The interface adapts the workstation data bus to RS-422.

The camera interface also supplies the CCD camera opto couplers with power.

Control of MAMMOMAT and biopsy unit

The workstation controls the MAMMOMAT via RS-232 communication.

Biopsy calculations are performed by the workstation and target coordinates are sent to the biopsy unit.

Monitor

NOTICE

If a 21" monitor is installed be sure not to change the proportions between length and width. If this has been done anyway, press the auto sizing and centring button marked with a symbol (Fig. 4) to restore the previous settings.

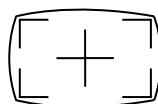


Fig. 4 Symbol

Power supply workstation main unit

NOTICE

When setting the main voltage on the contact, relate to local main voltage.

CCD camera

The CCD camera consists of a printed circuit board, CCD sandwich and RS-422 interface. When the x-ray beams reach the CCD sandwich, electric energy proportional to the x-ray energy is produced. The produced electrical current is A/D converted to digital information and sent to the workstation RAM via an RS-422 interface.

Biopsy unit

The biopsy unit is used for performing biopsy examinations.

The biopsy unit consists of a needle positioning device which is firmly attached to an 18 cm x 24 cm object table with a cut-out contour that is superimposed on the CCD. External stereo diaphragm, needle support and compression plate are also included in the biopsy unit. The biopsy unit can easily be attached to the swivel arm of the MAMMOMAT.

There are three displays on the biopsy unit front. The displays show spatial deviation of needle tip from the suspect point in x-, y- and z-axis calculated by the workstation.

Three potentiometers monitors the actual position of the needle tip.

This page intentionally left blank.

CAUTION

The CCD camera has to be handled with extreme care, it is very sensitive to mechanical shocks and temperature. When not connected do not touch the pins in the camera contacts. Shock and temperature sensors are integrated in the camera.

CAUTION

The CCD camera is sensitive to mechanical shock and shall always be stored in the attaché case, delivered with the system, when disconnected from the biopsy controller.

The camera shall be used within 10–30° C.

The camera shall be transported or stored within 0–40° C.

This page intentionally left blank.

It is very important that any intervention in the equipment shall start with disconnecting it from the power supply with the main circuit breaker. To prevent accidental triggering of high voltage and radiation, set the switch S2 (SS) on board D702 to OFF (lower position, no triggering of the SS relay).

CAUTION

Before switching off the power for the workstation it is important to first shut down the application software (see Logout and turning off the workstation in Supplement to the Instructions for Use MAMMOMAT 3000 - Opdima®).

⚠ WARNING

If the system is only switched off at the control panel or with S2/ D711 in the MAMMOMAT generator, line voltage will still be present at the generator line connection, line filter Z1, Z2, transformer T1, transformer T10 and board D711 (see MAMMOMAT 1000/3000 Nova Wiring Diagram). The Opdima® is switched off separately.

⚠ WARNING

After shut-down of the system, there may still be about 380 V DC present on the intermediate circuit of the MAMMOMAT generator. This will be indicated by LED V24 on board D710. The voltage will drop to less than 30 V within about 3 minutes, the LED goes out at about 30 V.

CAUTION

The boards contain electrostatic highly sensitive components requiring particular care in their handling (ground before making contact and place only on a conductive surface, see Technical Information TI 219). After completed hardware service, a protective earth measurement shall be performed.

This page intentionally left blank.

General

There are three different user levels in the Opdimas[®] software.

- Regular user
- Administration user
- Service user

Administration users use the Service mode to perform the following tasks:

- Calibration of the biopsy unit
- Setting up needles
- User administration
- Network setup, DICOM nodes
- Backup functions

Advanced service

Service users can additionally access advanced service functions. To get access you have to log in with the user name “service”. The password is obtained from Siemens Service Centre.

The following advanced service functions are included:

- Country settings
- Unit tests
- Software upgrade
- Restore disk
- Camera calibration/maintenance
- Disk cache
- Printer setup
- Log administration
- View log

NOTICE

The Opdimas[®] external diaphragm must be used for tests and calibrations involving radiation.

Selection of mode

After a successful login, the mode selection dialog is displayed.

1. Press the Service button.

Press Logout to cancel.

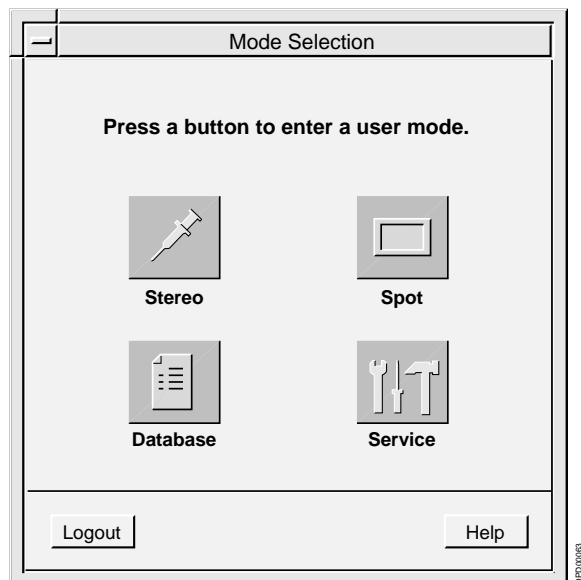


Fig. 1 Mode selection dialog

2. Select desired service function from the service dialog.

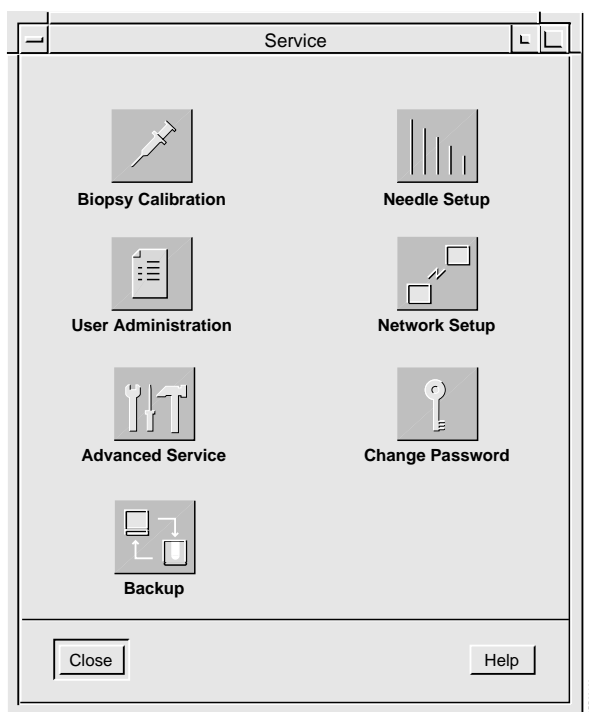


Fig. 2 Service dialog

Calibration of the biopsy unit

Calibration of the biopsy unit is carried out by using a stereo calibration phantom with targets at fixed positions. By releasing two stereo pair exposures with a fine needle adjusted to different targets, the system is automatically calibrated. Follow the instructions given on the screen during the calibration procedure.

1. Press the Biopsy Calibration button in the service dialog.
The following is shown on the screen:

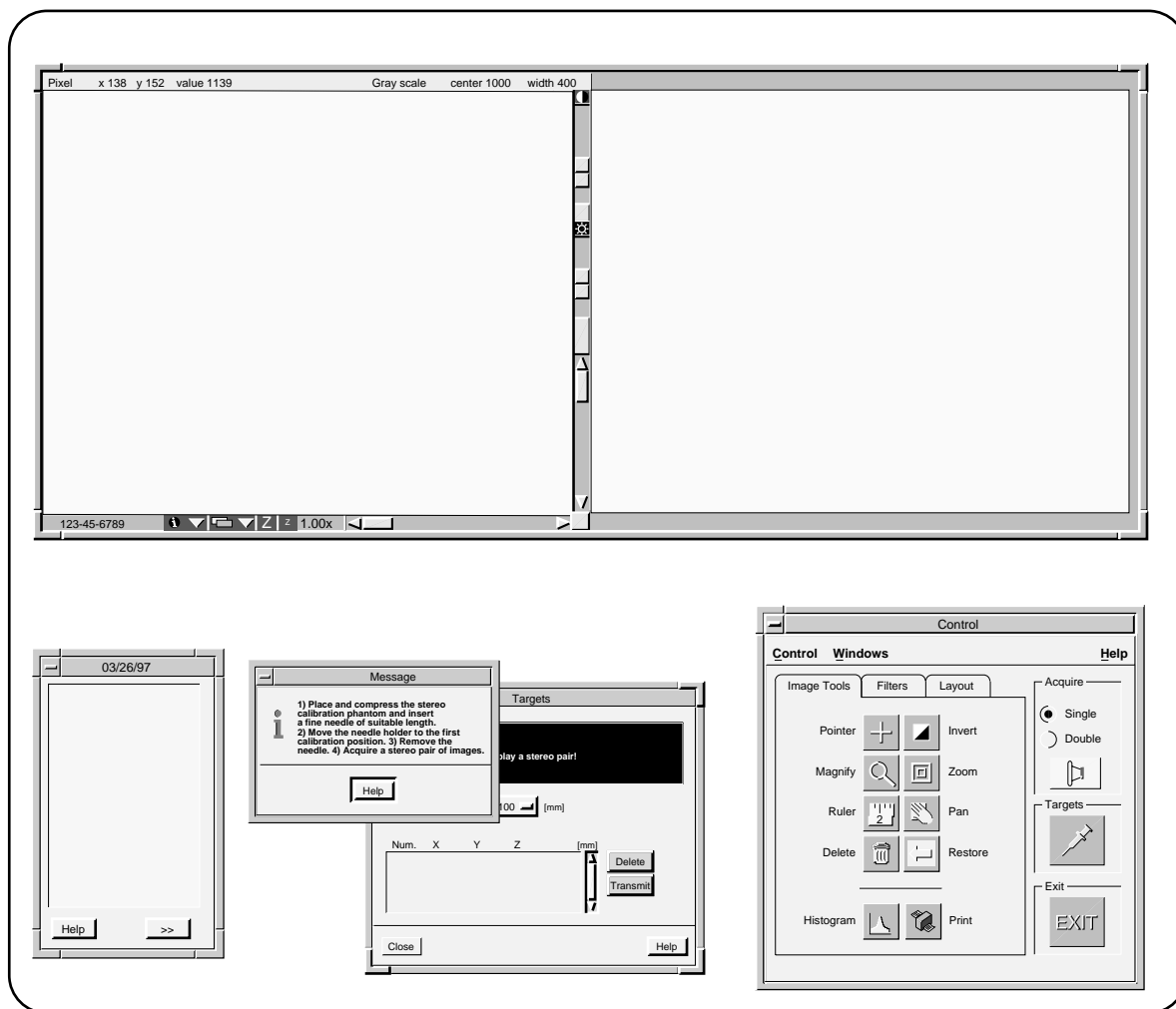


Fig. 3 Dialogs during calibration

2. Place the stereo calibration phantom on the biopsy table and compress.

With the targets facing the patient side, the phantom fits in the opening of the stereo compression plate, see Fig. 4. Target 2 shall be positioned 25–28 mm in negative x-axis direction and 12 mm from the object table side.

3. Select a fine needle to be used and measure the length.

Measuring the needle length is described on Page 5 - 7.

4. Choose Other in the Needle option menu and enter the selected needle length (minimum 90 mm).

OPD00385

5. Insert needle guides corresponding with the needle diameter.
6. Insert the fine needle into the needle guides.
7. Move the needle to Target 1 of the stereo calibration phantom, see Fig. 4, by using the adjustment knobs. Position the needle according to Fig. 5.

NOTICE

Make sure the needle is moved to Target 1.

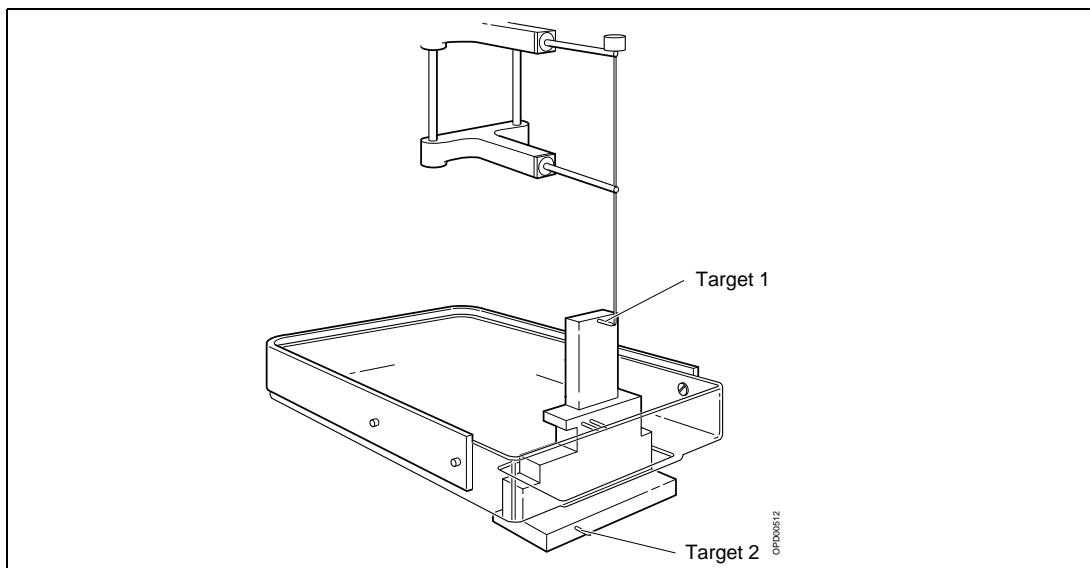


Fig. 4 Stereo calibration phantom

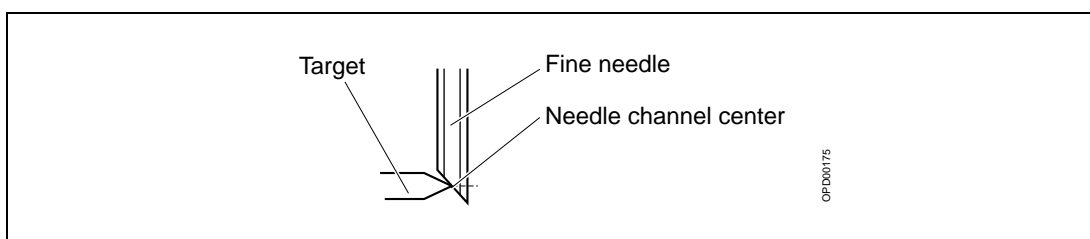


Fig. 5 Position of needle

8. Remove the needle.
9. Acquire and release a stereo pair of images.
Set the exposure parameters to 25 kV and 28 mAs in manual mode.
10. Check the reference marks and adjust if necessary, mark Target 1 and press Transmit.

NOTICE

When performing the biopsy calibration set the magnification to 0.7 and change Contrast/brightness in order to find all targets in the phantom. Be sure to mark the target at which the needle tip was positioned.

11. Insert the fine needle into the needle support of the biopsy unit.
12. Move the needle to Target 2 of the stereo calibration phantom by using the adjustment knobs.

13. Remove the needle.
14. Acquire and release a stereo pair of images.
Set the exposure parameters to 25 kV and 28 mAs in manual mode.
15. Check the reference marks and adjust if necessary, mark Target 2 and press Transmit.

When the calibration is successfully calibrated, the following message is displayed:

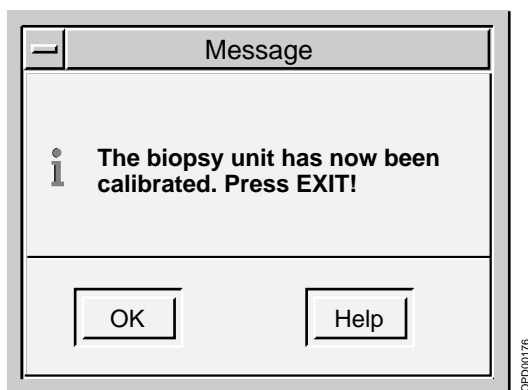


Fig. 6 Information message

16. Perform a final check according to Verifying the calibration of the biopsy unit on Page 10 - 1 to make sure that the biopsy unit works properly.

Setting up needles

1. Press the Needle Setup button in the service dialog.

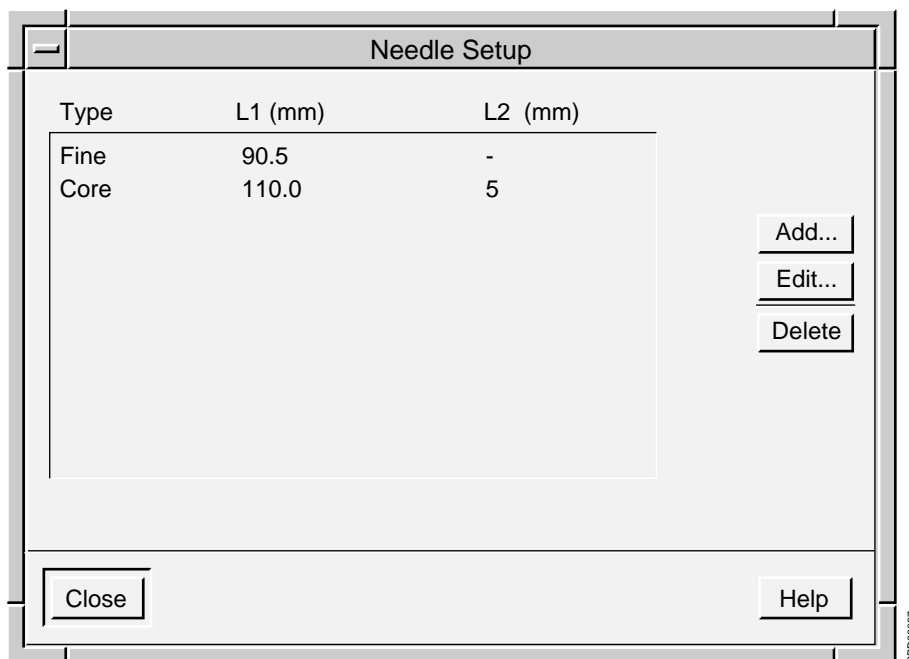


Fig. 7 Setting up needles dialog

2. Press Add... .

You can only add needles which match the needle guides with the fixed diameters: 0.7, 0.9, 1.2, 1.65, or 2.1 mm.

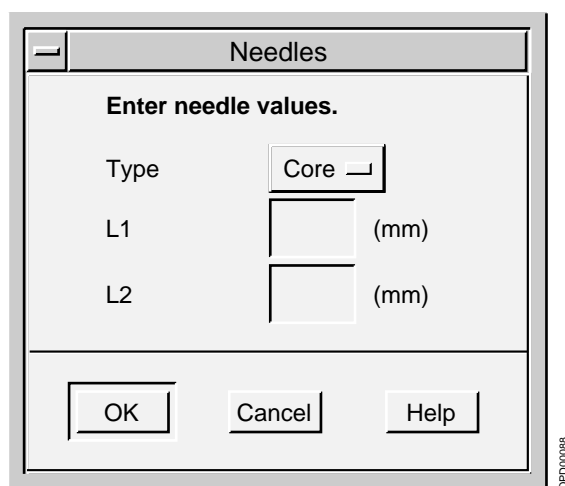


Fig. 8 Needle values dialog

3. Select needle type (fine or core) from the option menu and enter needle length in the text field.

L1 shall be 30 to 175 mm when selecting fine needle, L2 is not applicable for fine needle. When selecting core, L1 shall be set to <170 mm.

CAUTION

Ensure by measuring the length of the specified needle or core gun with needle that the right values are entered.

If the stroke length of the core gun is changed, a new needle has to be selected from the Targets dialog. The core needle length L1 shall be measured with the core needle in outer position and mounted in the core gun.

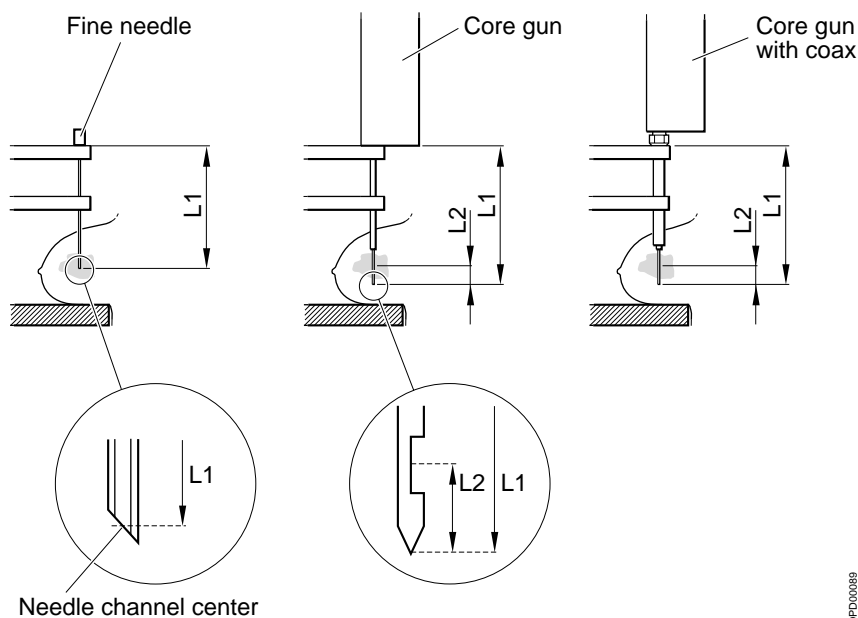


Fig. 9 Measuring of needle length

NOTICE

When using core gun, the system safety margin to avoid hitting the biopsy table is 5 mm.

Editing needle values

1. Select a needle from the needle setup dialog.
2. Press Edit... .

Backup of temporary storage media

As a service user it is possible to make a complete backup of all folders existing on the temporary storage media (the hard disk).

1. Press the Backup button in the service dialog.

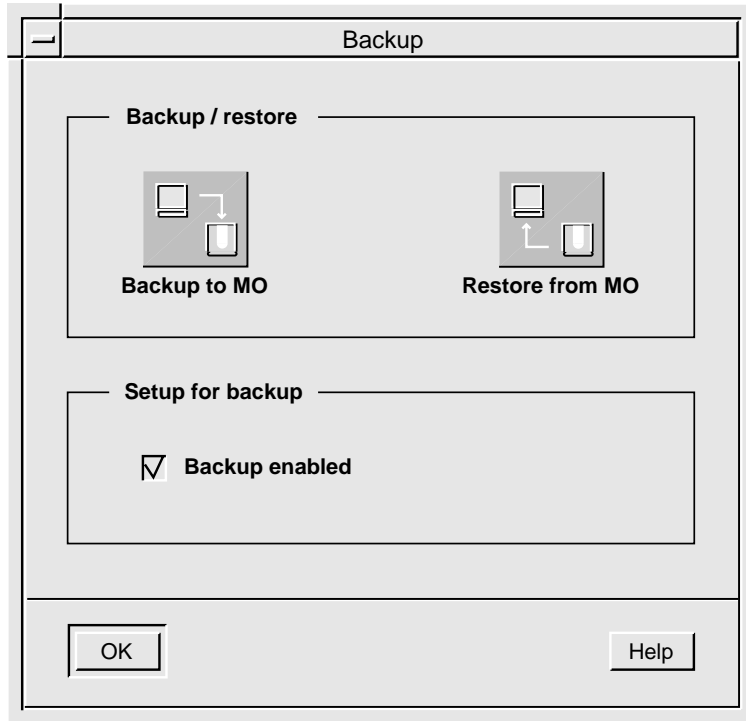


Fig. 10 Backup dialog

2. In the setup for backup section disable the backup.
3. Press the Backup to MO button.
4. Follow the instruction given in the message dialogs.
5. Enable the backup in the backup dialog.

WARNING

Backup of temporary storage media cannot replace the regular use of backup described in Using the Backup function in Supplement to the Instructions for Use MAMMOMAT 3000 - Opdima®.

Backup of temporary storage media only saves the content present on the hard disk at the actual moment. Images may have been erased by the disk cache system (see Page 5 - 23) or software reinstallation.

NOTICE

This backup might take a long time and require several MO disks since all folders on hard disk will be copied to backup MO.

Use of Advanced service functions

General

To get access to the advanced service menu you have to log in with the user name “service”. The password is obtained from Siemens Service Centre.

The advanced service is carried out by using a number of text dialogs. Selection of an item in a dialog can be done in three different ways:

- Use the up/down arrows of the keyboard to step through the fields (underlined)
- Use the space bar to step through the fields
- Type the corresponding number or letter of the item

Press Enter to execute.

Main menu

The main menu is used to access all the advanced service functions.

1. Press the Advanced service button in the service dialog.
2. Select a function from the menu.

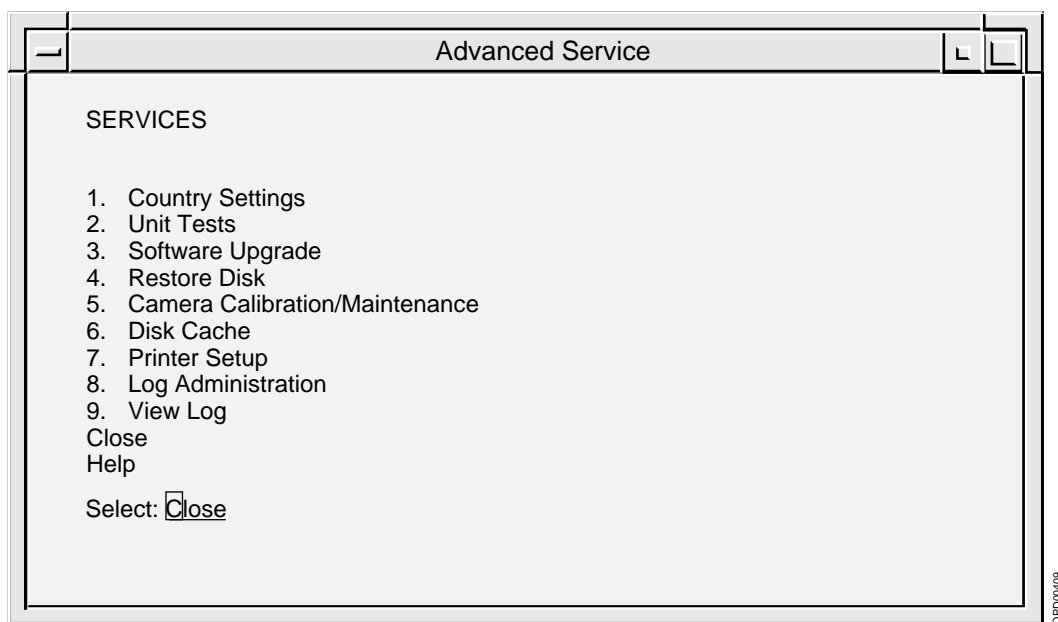


Fig. 11 Advanced service dialog

Country Settings	Function for selection of country of installation (your country).
Unit Tests	Function for performing tests of subassemblies.
Software Upgrade	Function for performing upgrade of software and for modifications of software.
Restore Disk	Function for restoring data from MO disk.
Camera Calibration/ Maintenance	Function for performing installation and maintenance of CCD camera.

Disk Cache	Function for disk cache settings.
Printer Setup	Function for setup of printer.
Log Administration	Function for log settings.
View Log	Function for inspection of logs.
Close	Function for leaving the advanced service mode.

Country settings

National parameters are set in the country settings dialog. On delivery it contains default settings.

Advanced Service

COUNTRY SETTINGS

Language: English

Time Zone: MET

Time Settings: Year: 2000 Month: 10 Day: 24 Hour: 09 Minute: 21

Id Pattern: 111-11-1111

Date Pattern: m/d/y

Service Center: Undefined

Institution Name: Undefined

Department Name: Undefined

Action: Close Close/Help/Apply

Fig. 12 Country settings dialog

Available Languages	Shows the languages currently available by the system.
Language	Selection of language in the dialogs. Select with space bar. If language is changed, log out and log in to make the change take effect.
Time Zone	A selection of the time zone used by the system. Select with space bar. For faster selection, type the first letter of the desired time zone name.
Time Settings	Value of system clock. If changed, Apply will reboot. NOTICE! The database is updated every time an examination is performed or when an image is loaded from a MO disk. The system clock cannot be set to a point earlier than the last update.
ID Pattern	Selection of ID-number structure. Type pattern with the keyboard: "1" for digit, "a" for letter and an arbitrary character for punctuation mark. For example: 111-11-1111.
Date Pattern	Selection of date structure. Type pattern with the keyboard: "y" for year, "m" for month, "d" for day and an arbitrary character for punctuation mark. For example: m/d/y or y-m-d (do not type mm/dd/yy or yy-mm-dd).
Service Center	Type the appropriate service center.
Institution Name	Type the appropriate name for the institution.
Department Name	Type the appropriate name for the department.

1. Select Country Settings in the advanced service dialog.
2. Set the values for Language, Time Zone, ID Pattern, Date Pattern, Service Center, Institution Name and Department Name in each respective fields.
3. Select Apply and press Enter to execute the changes.

Test of units

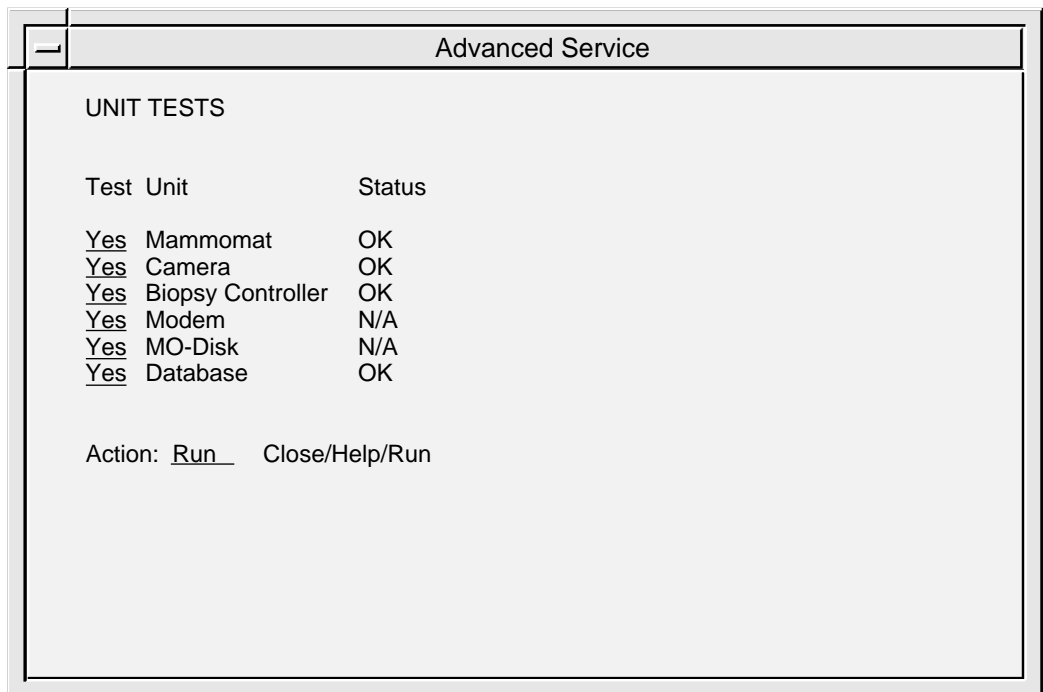


Fig. 13 Unit tests dialog

NOTICE

Before performing unit tests, the biopsy unit needs to be mounted to the MAMMOMAT. Otherwise the biopsy controller test will fail.

1. Select Unit Tests in the advanced service dialog.
2. Select Yes or No to select/deselect the parts to be included in the test.
3. Select Run and press Enter.

Each unit that passes the test will be indicated with an OK message. If it does not pass, a Failed message will appear. Units that have not been tested are indicated with an Untested message.

The units are tested according to the following:

- MAMMOMAT - test if there is a connection and if the power is on
- CCD camera - test if there is a connection and if the power is on
- Biopsy controller - test if there is a connection and if the power is on
- Database - reading and writing in the database

Software upgrade

The software upgrade dialog is used to upgrade the current software version with a new version from CD-ROM and to do modifications of the software.

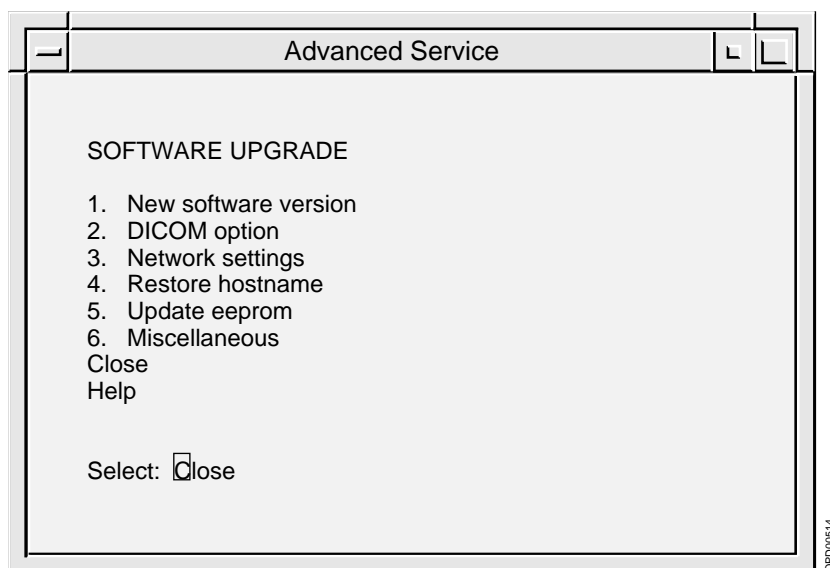


Fig. 14 Software upgrade dialog

New software version	Upgrade with new version from CD-ROM.
DICOM option	Enabling the DICOM option. NOTICE! A license key is necessary to install DICOM.
Network settings	Function for defining the network settings for the Opdimas [®] system.
Restore hostname	The hostname of the original workstation can be restored to a new workstation. Restore hostname from most recently used local/backup MO disk.
Update eeprom	For Ultra 1 only: Save system specific information to NVRAM.
Miscellaneous	Database and MO disk utilities.

New software version**NOTICE**

This menu will be used when installing the next version of software.

CAUTION

Make sure that all examinations are stored on the MO disk before upgrading. The examinations are stored when exiting an examination or database session, see Storing data on MO disk in Supplement to the Instructions for Use MAMMOMAT 3000 - Opdima®.

1. Select Advanced service in the Service mode.
2. Insert the CD-ROM for the new ASW software version.
3. Select Software upgrade and New software version.



Fig. 15 New software version, advanced service dialog

4. Select Install.

Enabling the DICOM option

To enable the DICOM option you need the DICOM license key that is received when purchasing the DICOM option.

The DICOM licence is written on the invoice and the dispatch note. If not, please contact Siemens-Eléma AB mammography logistics department in Solna, Sweden for further information. The mammography logistics department will ask for serial No. and host ID of the workstation in question.

Instruction for obtaining the systems host ID and hostname

1. Switch on the workstation.
2. Log in as service user.
3. Press the **Control, Alt, Shift** and **!** keys at the keyboard simultaneously.
4. Switch to default behavior?
Press OK. The screen will flash for a second.
5. Place the mouse pointer on the screen background and press the rightmost mouse button.
6. Select New window.
7. Type "hostname" and press Enter.

NOTICE

Note the hostname (opdxxxx) shown on the screen.

8. Type "hostid" and press Enter.

NOTICE

Note the host id (xxxxxxxx) shown on the screen.

9. Type "exit" and press Enter.
10. Press the **Control, Alt, Shift** and **!** keys at the keyboard simultaneously.
11. Switch to custom behavior?
Press OK. The screen will flash for a second.
12. Done!

DICOM option

1. Select Advanced service in the Service mode.
2. Select Software Upgrade.
3. Select DICOM option.

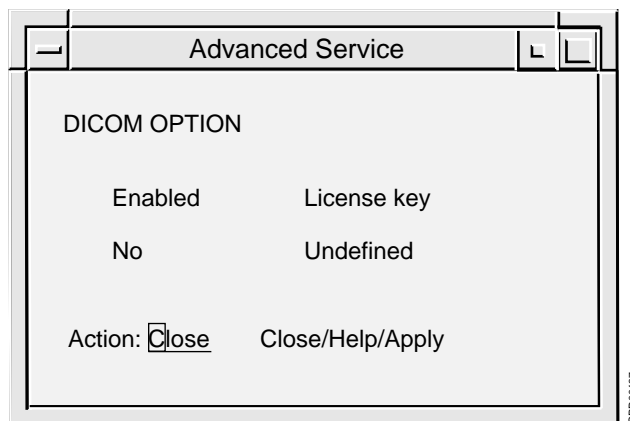


Fig. 16 DICOM option, advanced service dialog

4. Select Yes.
5. Enter License key and Apply.
6. Select y to reboot.
7. Configure DICOM nodes according to Network Setup in Supplement to the Instructions for Use MAMMOMAT 3000 - Opdima®.

Network settings

Enable network

1. Select Advanced service in the Service mode.
2. Select Software upgrade and Network settings.
3. Supply appropriate values for the network settings into the fields and apply.
(The values shall be supplied by the network administrator at the hospital.)

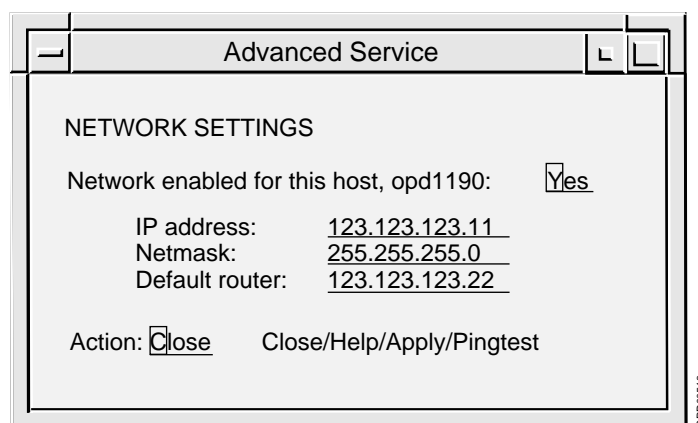


Fig. 17 Network settings

IP address	Address for the system in the network, written as four decimal numbers separated by periods, e.g. 123.123.123.11. No initial 0, e.g. 10.10.10.10 is OK but 010.010.010.01 is not OK.
Netmask	Netmask used in the network if IP standard subnetting is used.
Default router	The IP address of the default router, if this is used in the network. No initial 0, e.g. 10.10.10.10 is OK but 010.010.010.01 is not OK.

4. Select y to enable the network. The workstation will be turned off.
5. Connect the network cable.
6. Switch on the workstation.

Pingtest

A pingtest checks that the default router is alive, it can be performed when the network is enabled.

1. Select Pingtest in the Network settings dialog.

NOTICE

Some alive routers do not reply to pingtest, check with network administrator.

To test if a host is alive or not, write the IP address temporarily in the field for Default router and perform a pingtest.

Miscellaneous

1. Select Advanced service in the Service mode.
2. Select Miscellaneous.

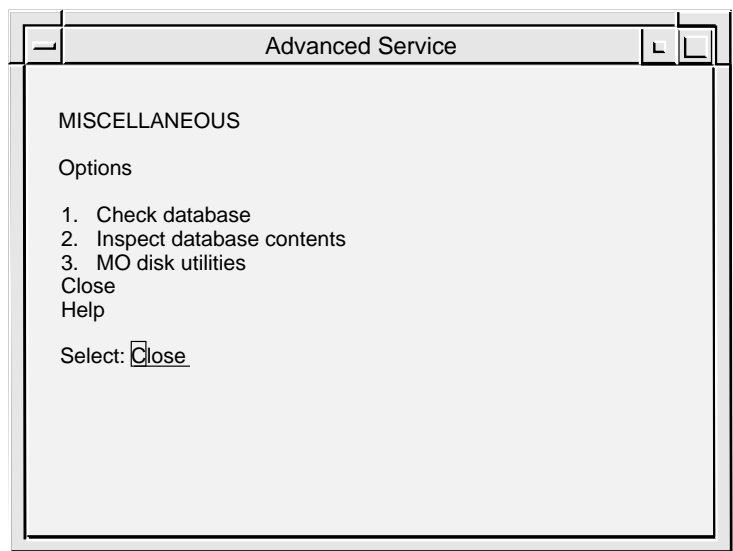


Fig. 18 Miscellaneous, advanced service dialog

Check database

This will check and adjust invalid folder names that may have been introduced in earlier software versions.

1. Select check database.

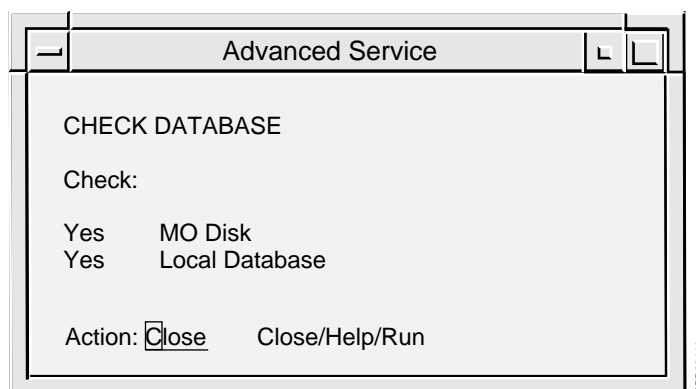


Fig. 19 Check database, advanced service dialog

2. Select the database to check, MO Disk and/or Local Database.
3. Select Run.

Inspect database contents

1. Select Inspect database contents.
A Text Editor window with a log file will appear on the screen, see example in Appendix 2.

NOTICE

**The log file is only a print-out of the data in the database.
Changes in the log file will NOT affect the database.**

2. When done, close the Text Editor window.


MO disk utilities

Fig. 20 MO disk utilities, advanced service dialog

Check MO disk	Check MO disk file system integrity.
Repair MO disk	Repair MO disk file system. NOTICE! Run only if Check MO disk shows that something has to be done.
View output	Lists output from latest "Repair MO disk" execution.
Identify MO disk	Shows the kind of disk (local storage/backup/transfer...).
List MO directories and files	View contents of MO disk.
Create a clean MO disk	Format a new disk or erase contents of a formerly used disk. WARNING! Use with extreme caution. Erases all content on MO disk.

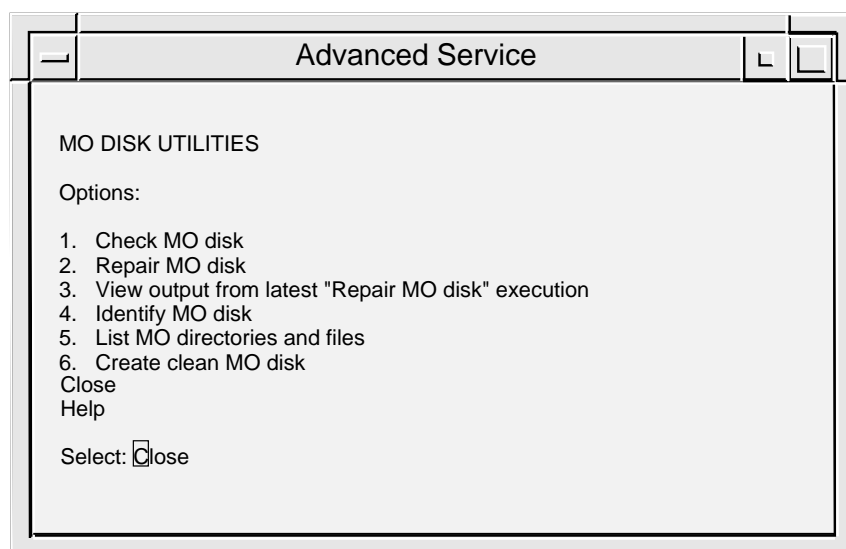Restoring data from MO disk

Copies the database from the MO disk to the hard disk. It is possible to use a local storage MO disk or a backup MO disk.

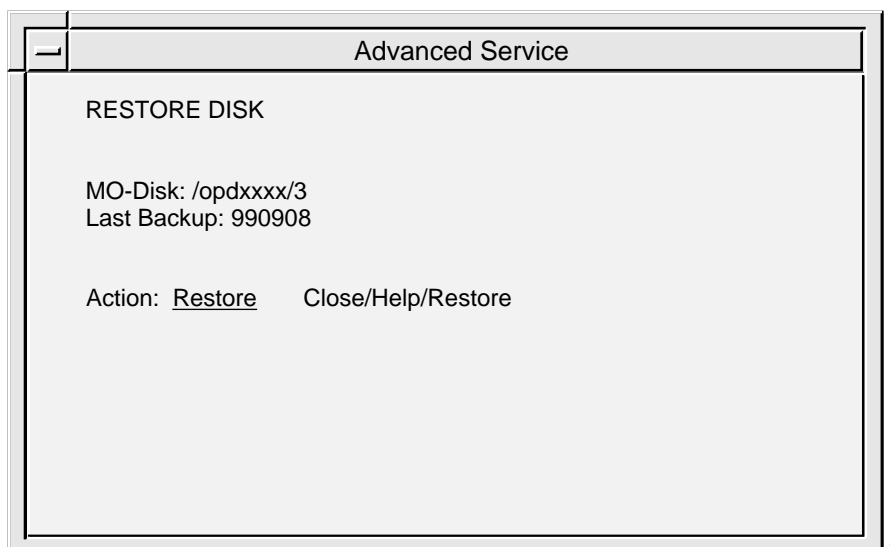


Fig. 21 Restore disk dialog

MO disk	Name of the inserted MO disk.
Last backup	Last date for storing on the MO disk.

1. Insert the most recently used MO disk in the MO unit.
2. Select Restore Disk in the advanced service dialog.
3. Select Restore and press Enter to copy the data to the hard disk.



WARNING

When restoring data, be sure to insert the MO disk that was used most recently before the reinstallation of software.

Date for last update of MO disk is displayed on the monitor before confirmation of restore.

If the system is restored from an older MO disk, the most recent folders will be lost from the database and the new MO disk numbering can be incorrect. Please contact HSC for more information.

NOTICE

Set up the printer once again after software reinstallation. Printer info can not be restored.

CCD camera calibration and maintenance

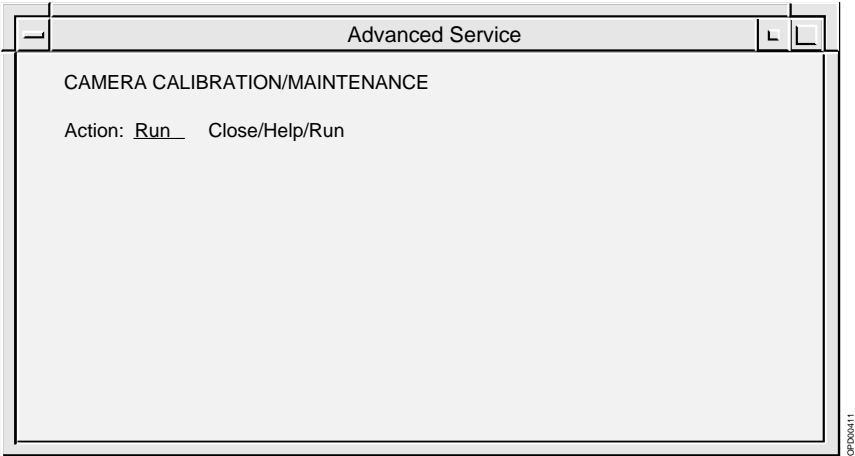


Fig. 22 Camera calibration/maintenance

NOTICE

This calibration shall be performed with the object table (grid or non grid) that the customer is going to use in most cases.

NOTICE

Do not use the stereo table when calibrating.

Fill in your values in the test protocol for CCD camera calibration (see Appendix 3).

Grid table

- 1. To avoid grid lines in the calibration images, temporary increase the grid speed during calibration as follows:
Change the grid fast speed time to 1500 ms, the grid fast speed to 99% of max and the grid slow speed to 40% of max (Fig. 23) using the service PC.

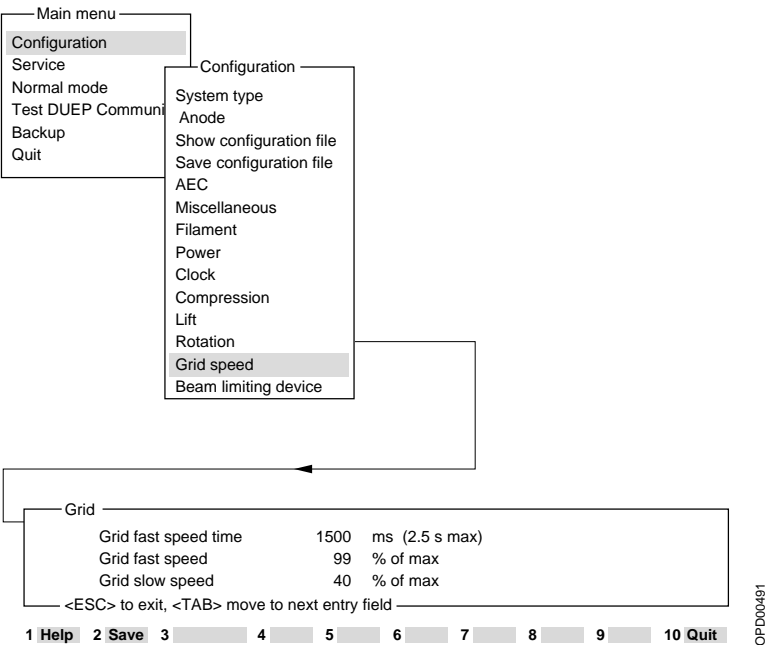


Fig. 23

2. Place 40 mm acrylic plastic (PMMA) on the table, covering the CCD.
3. Select Camera Calibration/Maintenance in the advanced service dialog.
4. Select Run and set 27 kV, Mo/Mo and AEC mode, on the generator.
Follow the instructions displayed on the screen.

Between each exposure the acrylic plastic shall be moved slightly (the reason being that impurities in the plastic shall not influence the calibration).
After the exposure series, correction tables will be calculated. This will take 10 minutes at most.

After a successful calibration (10 exposures in each mode - normal resolution and high) the message "Camera successfully calibrated" will appear on the screen.

5. Change the grid speed back to original values.
6. Perform check of Opdima® AEC and check of resolution. See Check of Opdima AEC on Page 10 - 2 and Check of resolution on Page 10 - 4.

Non grid table

1. Place 30 mm acrylic plastic (PMMA) on the table, covering the CCD.
2. Select Camera Calibration/Maintenance in the advanced service dialog.
3. Select Run and set 26 kV, Mo/Mo and AEC mode, on the generator.
Follow the instructions displayed on the screen.

Between each exposure the acrylic plastic shall be moved slightly (the reason being that impurities in the plastic shall not influence the calibration).
After the exposure series, correction tables will be calculated. This will take 10 minutes at most.

After a successful calibration (10 exposures in each mode - normal resolution and high) the message "Camera successfully calibrated" will appear on the screen.

4. Perform check of Opdima® AEC and check of resolution. See Check of Opdima AEC on Page 10 - 2 and Check of resolution on Page 10 - 4.

Disk cache settings

Storage parameters can be set in the disk cache dialog:

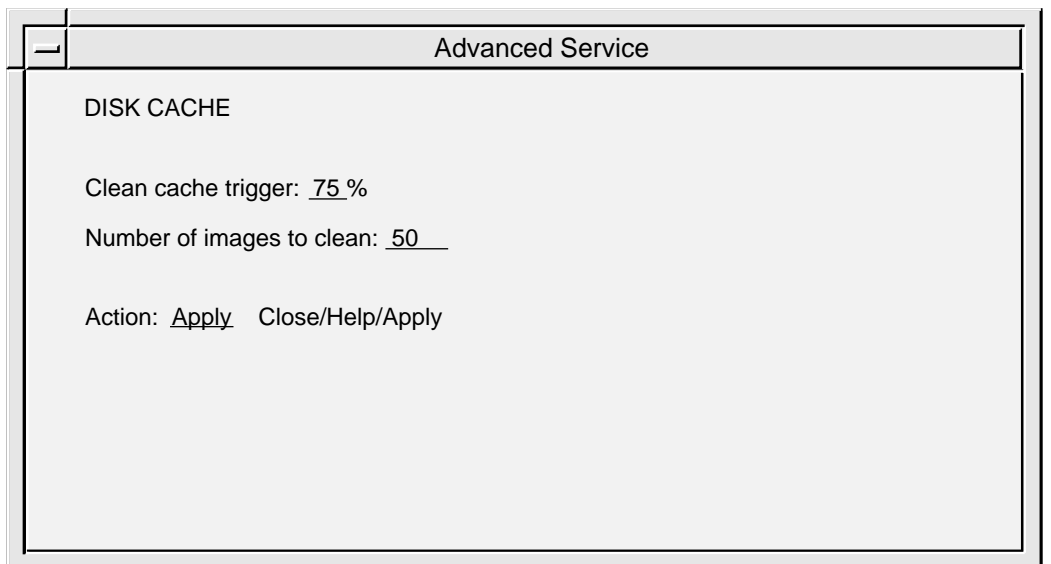


Fig. 24 Disk cache dialog

Clean cache trigger	Cleans a number of images from the cache memory (the temporary storage media i.e. hard disk) when it has been filled up to a percentage of the capacity. Normally set to 75%.
Number of images to clean	The number of images to erase once the trigger level has been reached. Normally set to 50.

1. Select Disk Cache in the advanced service dialog.
2. Set the values.
3. Select Apply and press Enter.

Every time an examination is performed a folder containing the image file, the image icons and the image header file is created. When an examination is finished the complete folder is copied to the MO disk. If a folder has been copied to a MO disk, its images will be erased from the hard disk if the cache limit is reached.

Use of printer setup

Stand alone Opdima system

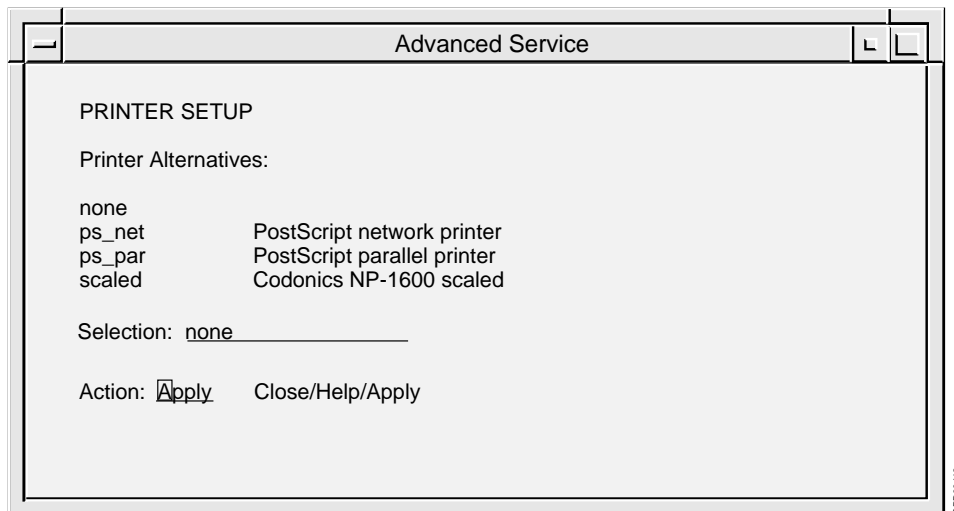
NOTICE
Connection of stand alone printer.


Fig. 25 Printer setup dialog

1. Select printer by pressing the Space button.
 - 1.1 ps_net: PostScript network printer, for postscript printer connected to the network connector of the workstation see 7/Fig. 10 on Page 6 - 11. The printer IP address has to be set to 10.10.10.2.
 - 1.2 ps_par: PostScript parallel printer, for postscript printer connected to the parallel port of the workstation see 8/Fig. 10. on Page 6 - 11.
 - 1.3 scaled: Codonics NP-1660 or NP-1600 scaled, for a Codonics printer connected to the network connector of the workstation see 7/Fig. 10 on Page 6 - 11. The printer IP address has to be set to 10.10.10.2.
2. Select Apply and press Enter.
3. Connect your selected printer after turning off the system.

NOTICE
If no printer is connected select none, apply and restart the system.

Networked Opdima system

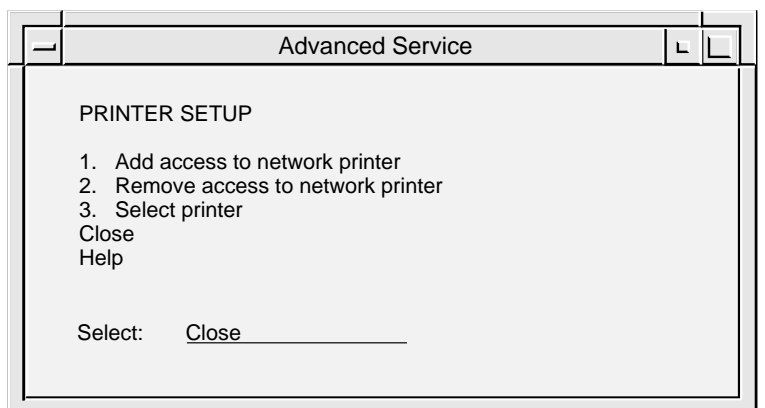
NOTICE**Connection of networked printer.**

Fig. 26 Printer setup dialog when connected to a network

When the Opdima® system is connected to a network (see Page 5 - 17) it can be configured to print on a printer connected to the network. Supported printer types are postscript printers and printers that can handle XWD files (X Window Dump).

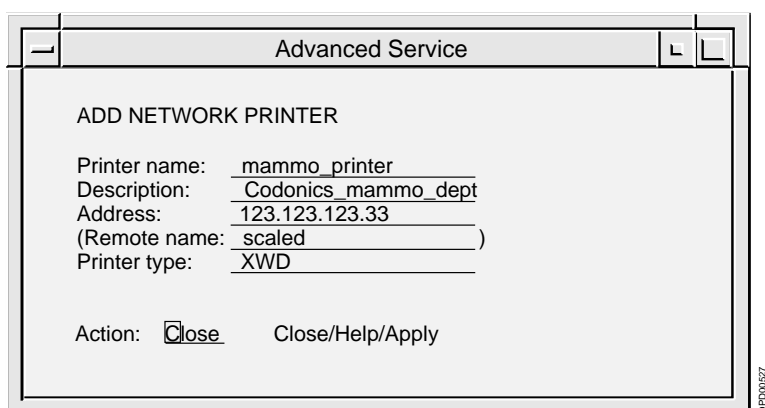
Add access to network printer

Fig. 27 Add network printer dialog

- Printer name: Logical name presented when selecting printer in Advanced service. No space characters are allowed. E.g. mammo_printer.
- Description: Textual description, for example the location of the printer. E.g. Codonics_mammo_dept.
- Address: IP address for printer, e.g. 123.123.123.33.
No initial 0, e.g. 10.10.10.10 is OK but 010.010.010.01 is not OK.

Remote name: Name of printer at remote system. No space characters are allowed.
NOTICE! The remote name “scaled” must be used when using a Codonics printer.

Printer type: PS or XWD. Use PS for all postscript printers.
NOTICE! XWD shall be used with Codonics printers.

Select printer

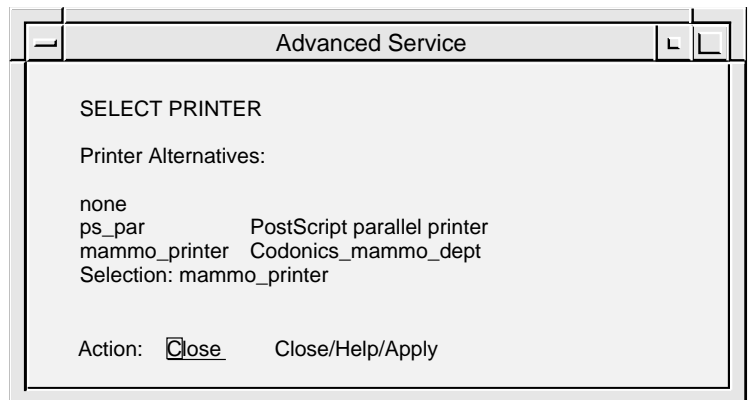


Fig. 28 Select printer dialog, example

Note that only ps-par and network printers that have been added can be selected when the network has been enabled for the system.

Remove access to network printer

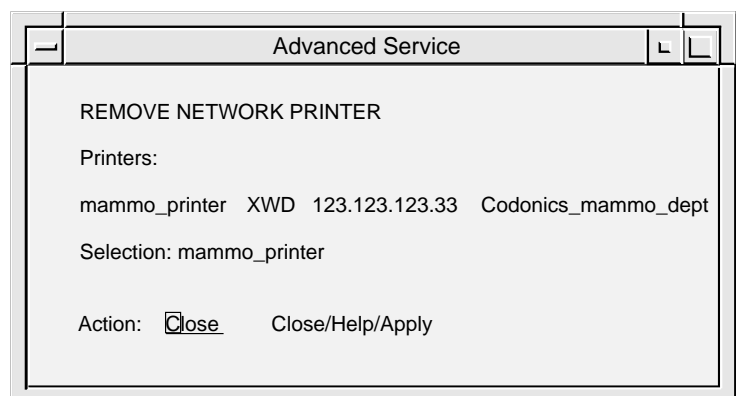


Fig. 29 Remove network printer dialog, example

Previously defined printers can be removed from the list of printers. Note that the printer that is currently selected cannot be removed.

Log administration

Log files are stored on the temporary storage media and on the MO disks. In the log administration dialog, the logging can be turned on/off.

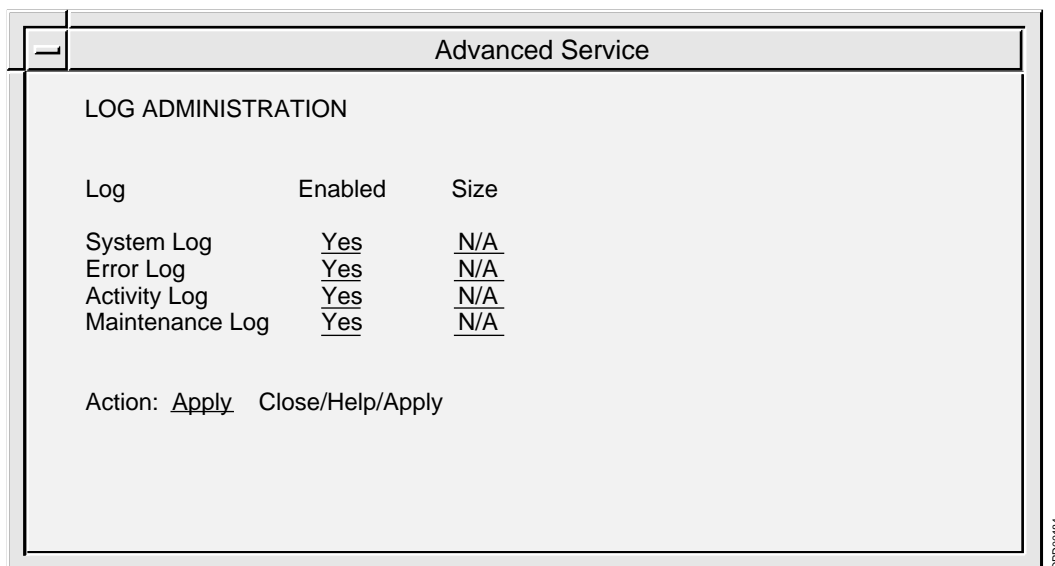


Fig. 30 Log administration dialog

1. Select Log Administration in the advanced service dialog.
2. Select Yes/No to turn on/off each respective log.
3. Select Apply and press Enter.

Log inspection

Through the view log dialog, the log files can be inspected.

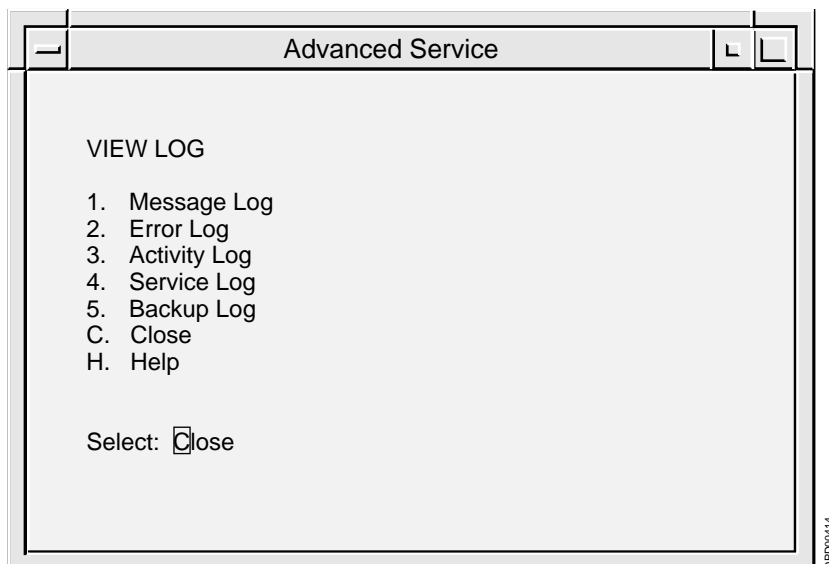


Fig. 31 View log dialog

Message Log	Log for messages presented to the user.
Error Log	Log for Opdima® related errors.
Activity Log	Log for all user activities.
Service Log	Log for activities in service mode.
Backup Log	Log for activities made in backup mode.
Close	Function for leaving the view log dialog.

1. Select View Log in the advanced service dialog.
2. Select a log from the menu in the view log dialog.

General

After performing hardware service, protective earth measurement has to be performed, see Protective earth measurement on Page 10 - 5.

WARNING

Switch off the power to the system at the main circuit breaker after complete power off.

CAUTION

If the workstation shall be powered off, use the “poweroff” command and make sure that the system has completed the process before switching off any main power. If power to the workstation is disabled during the power off sequence or without using the “poweroff” command, data could be lost.

For Ultra 10 only: After the “poweroff” command has been started the OK prompt will appear. To switch off the power completely press Control + power on key simultaneously.

MO unit (no CD)

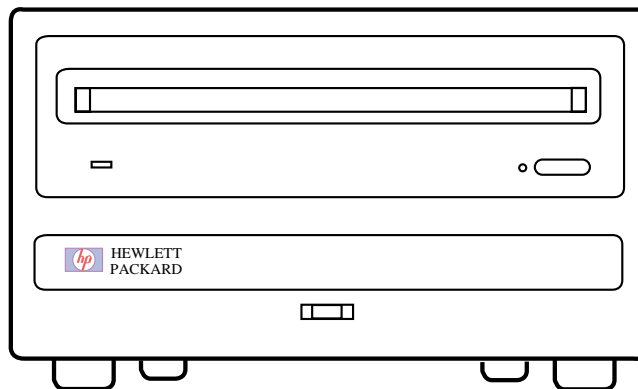


Fig. 1 The MO unit, front view

Removing the MO unit

NOTICE

If MO unit is defective, the complete unit must be replaced.

1. Disconnect the power supply cable (AC in) and SCSI cable from the back of the unit (Fig. 2).
2. Unplug the SCSI terminator from the back of the unit (Fig. 2).
3. Remove the unit.

Installing the MO unit

NOTICE

Make sure that the SCSI cables are properly connected and fixed to the units and that the cable is not bent to much just behind the MO unit. Also, check that the MO unit is terminated correctly. If not, this might generate a lot of problems e.g. failure to write images to MO disk, corrupted MO disk or corrupted hard disk.

SCSI ID	Disk drive type	Disk drive function
1	Magneto-optical drive	Permanent storage and backup media (MO disk)

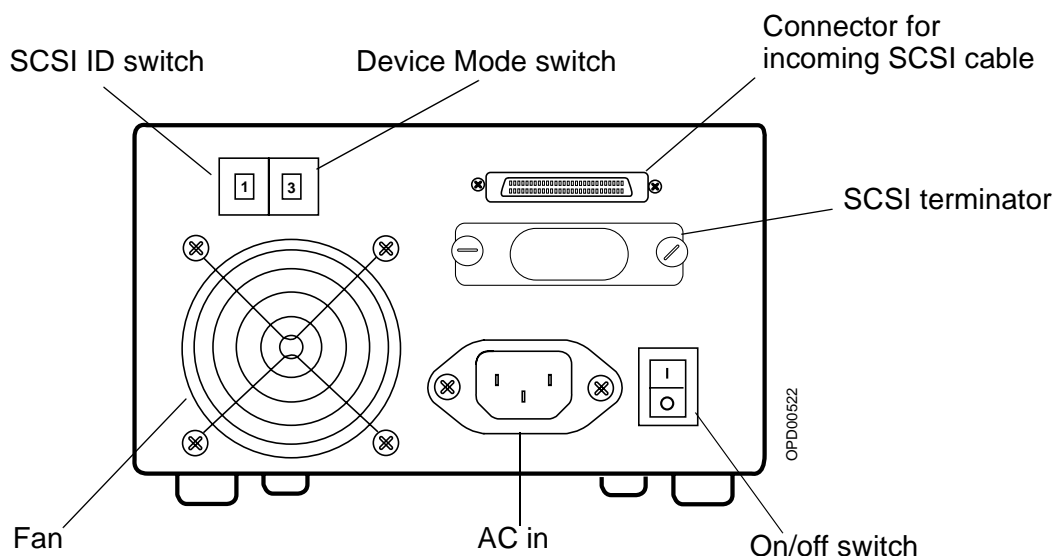


Fig. 2 MO unit, back view

1. Set SCSI ID switch to 1 (SCSI address).
2. Set Device mode switch to 3 (Direct access/verify on).
3. Connect the power supply cable, the SCSI terminator and the SCSI cable to the unit.

Tests and Adjustments

1. Switch on the system.
2. Ensure that the system is operating properly.

CD/MO unit

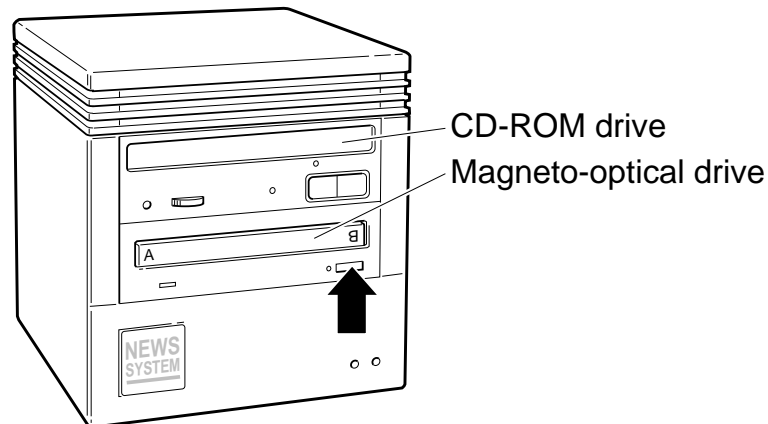


Fig. 3 CD/MO unit

Removing the CD/MO unit

NOTICE

If the CD drive or MO unit is defective, the complete unit must be replaced.

1. Disconnect the power supply cable (AC in) and SCSI cable from the back of the unit (Fig. 4).
2. Unplug the SCSI terminator from the back of the unit (Fig. 4).
3. Remove the unit.

Installing the CD/MO

NOTICE

Make sure that the SCSI cables are properly connected and fixed to the units and that the cable is not bent to much just behind the MO unit. Also, check that the MO unit is terminated correctly. If not, this might generate a lot of problems e.g. failure to write images to MO disk, corrupted MO disk or corrupted hard disk.

SCSI ID	Disk drive type	Disk drive function
1	Magneto-optical drive	Permanent storage and backup media (MO disk)
6	CD-ROM disk drive	Software distribution

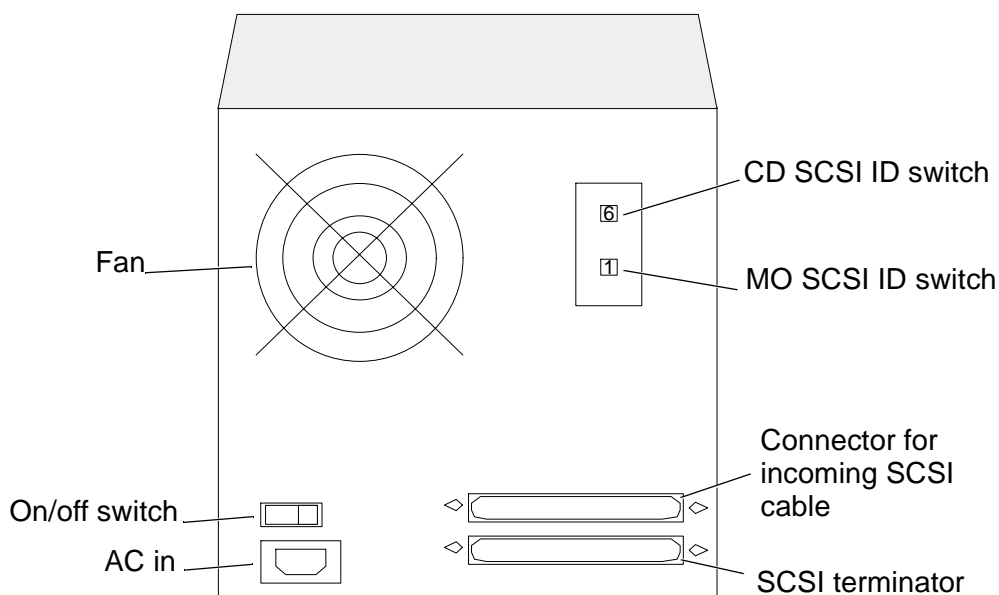


Fig. 4 CD/MO unit, back view

NOTICE

If the MO unit has no switches, the SCSI IDs are set internally.

1. Set CD SCSI ID switch to 6.
2. Set MO SCSI ID switch to 1.
3. Connect the power supply cable, the SCSI terminator and the SCSI cable to the unit.

Tests and Adjustments

1. Switch on the system.
2. Ensure that the system is operating properly.

Biopsy controller

Removal of biopsy controller

The biopsy controller is attached to the stand rear side. To remove the unit proceed as follows:

1. Remove the CCD camera and store it in its attaché case.

CAUTION

Do not touch the pins in the contacts of the camera.

The CCD camera is sensitive to mechanical shock and shall always be stored in the attaché case when disconnected from the biopsy controller.

2. Remove the biopsy controller cable duct.
3. Disconnect all cables.
4. Lift off the biopsy controller (1/Fig. 5) from the holders of the stand (2/Fig. 5).

WARNING

If the MAMMOMAT power is switched on, there is mains voltage present at the biopsy controller mains cord connector.

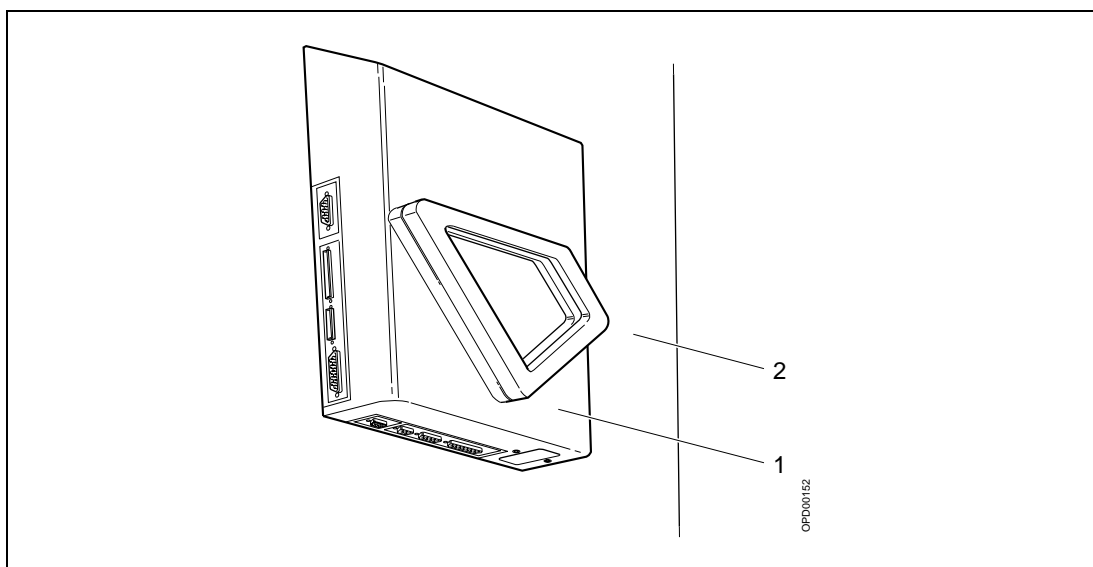


Fig. 5 Biopsy controller

Removal of biopsy controller cover

To remove the cover of the biopsy controller, proceed as follows:

1. Remove the 13 screws (4/Fig. 6) on the back side of the biopsy controller.
2. Remove the seven countersunk screws (2/Fig. 6) on the sides of the biopsy controller.
3. Lift off the cover, starting at the top.

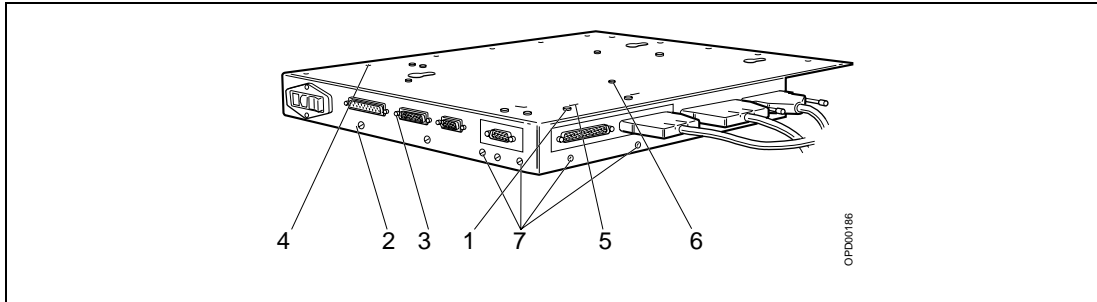


Fig. 6 Removal of cover

Replacement of components in the biopsy controller

Replacement of board D200

To replace the board D200 in the biopsy controller, proceed as follows:

1. Remove the screws (3/Fig. 6) holding the connectors X201, X202, X206, X207, X208.
2. Loosen the four screws (1/Fig. 6) on the back side of the biopsy controller.
3. Loosen the four screws (7/Fig. 6) on the side of the biopsy controller.
4. Remove the screws (6/Fig. 6) holding the board.
5. Remove the board and install the new board. If necessary, loosen the transformer.

Make sure to replace the isolation film at the connectors X203, X204 and X205.

When installing the board, pry the board holders carefully in position by using a screwdriver at the slots (5/Fig. 6) on the back side of the biopsy controller.

6. With an ohm meter verify that the isolation film at connectors X203, X204 and X205 is functioning properly. The resistance shall be in excess of 10 MΩ, when measuring between the back side of the biopsy controller and the shielded housing surrounding the pins of the connectors X203, X204 and X205. If the resistance is less than 10 MΩ, the isolation film has to be replaced.

NOTICE

Led V2 NVM error “Biopsy unit not responding” will appear after replacement of D200, see Biopsy unit not responding on Page 9 - 2.

Replacement of fuse in mains filter

The mains filter (3/Fig. 7) contains a fuse. To change fuse, proceed as follows:

1. Remove the fuse holder (1/Fig. 7) from the filter using a screwdriver and remove the old fuse (2/Fig. 7).
2. Install a new fuse in the fuse holder.
3. Insert the fuse holder into the filter according to Fig. 7.

CAUTION

The new fuse must be placed in the same position as the old one and the marks on the fuse holder and the mains filter must coincide, see Fig. 7.

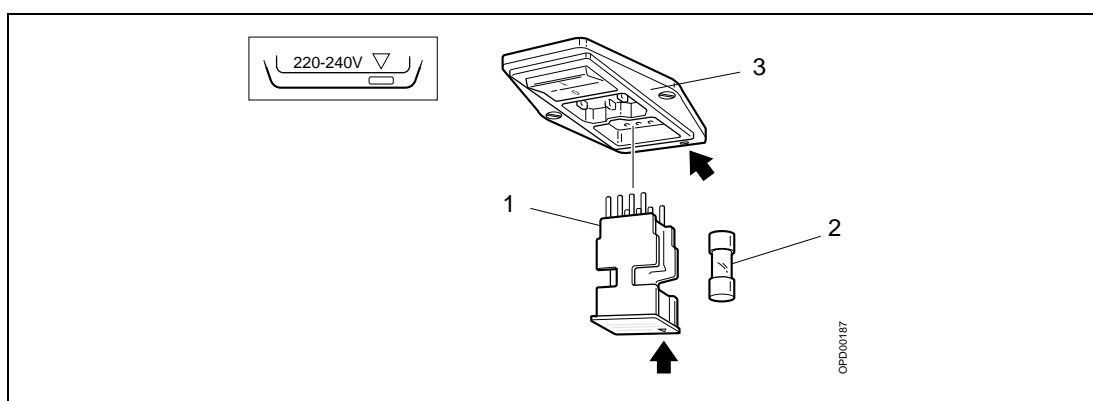


Fig. 7 Fuse in mains filter

Replacement of fuses on printed circuit board, D200

Fuses on the board D200 in the biopsy controller can be exchanged. The location of the fuses is shown in MAMMOMAT 1000/3000 Nova - Opdima® Wiring Diagram.

Replacement of EPROMs

To replace the EPROMs of the board D200 in the biopsy unit, proceed as follows:

1. Use PROM extractor to remove the EPROM.
2. Install the new EPROM.

Ultra 10 Cables

The biopsy controller is the central unit for the cable connections. The cables are connected to the biopsy controller according to the figure below:

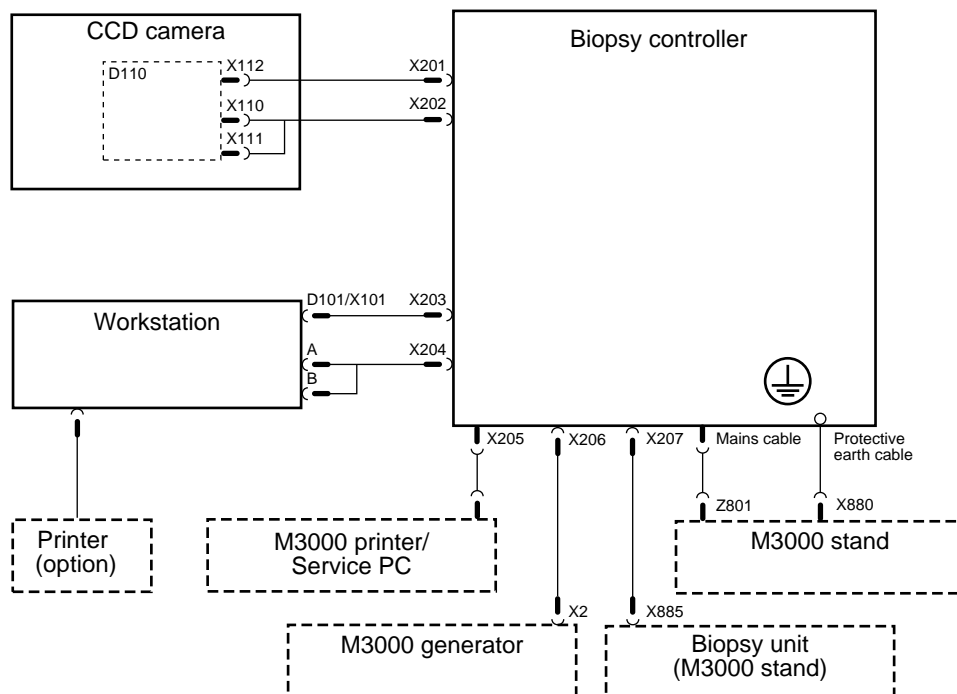


Fig. 8 Connection of cables

OPD00464

Ultra 1 Cables

The biopsy controller is the central unit for the cable connections. The cables are connected to the biopsy controller according to the figure below:

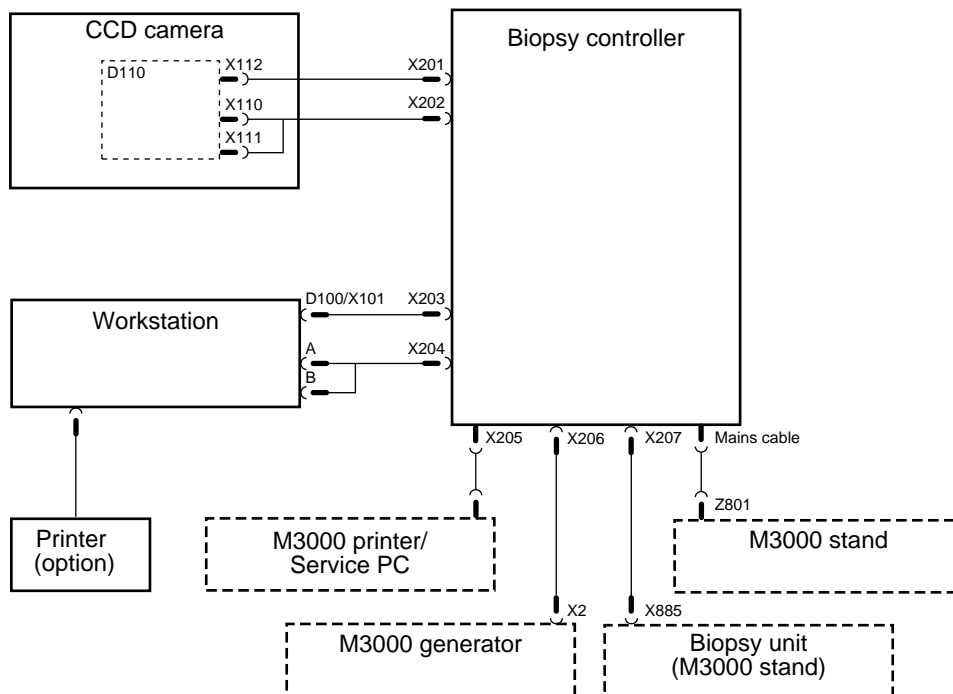


Fig. 9 Connection of cables

OPD00188

Biopsy unit

After any repair on the biopsy unit the accuracy has to be tested with the stereo calibration phantom and if necessary a calibration of the biopsy unit has to be performed.

Ultra 10 Workstation

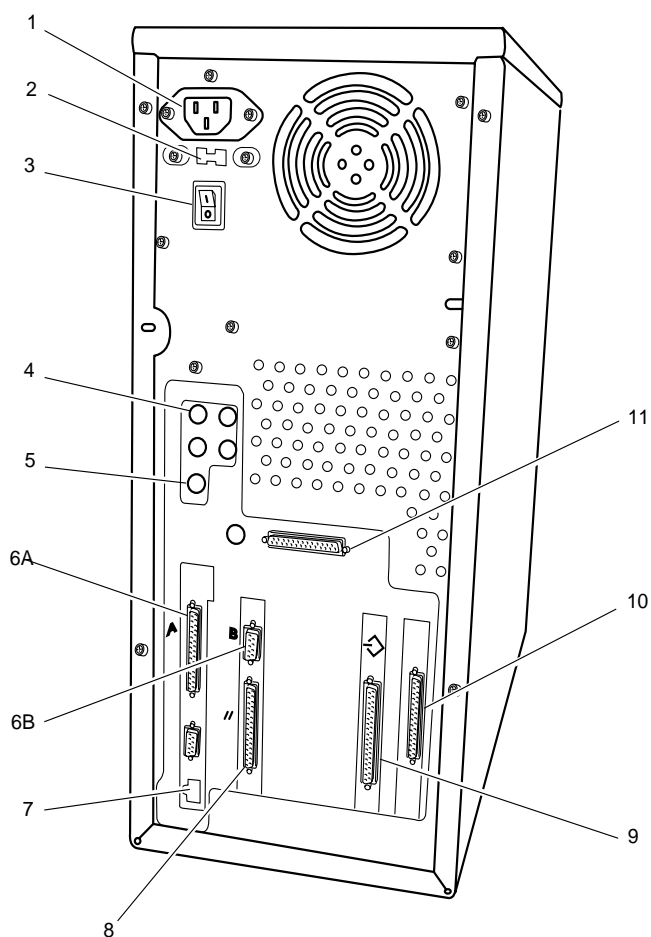
Connectors

CAUTION

Connecting the cable to the wrong connector may cause severe damage to the equipment.

CAUTION

If the workstation table (option) is used, the main unit shall be placed upright in the designated shelf on the framework of the workstation table (option).



- | | |
|--|--|
| 1 Mains power supply | 7 Network connector |
| 2 Line voltage selector | 8 Parallel port to local printer |
| 3 Power on / off | 9 SCSI interface |
| 4 MIC | 10 Camera interface via biopsy controller (X203) |
| 5 Keyboard interface | 11 Monitor interface |
| 6A Mammomat interface via biopsy controller (X204) | |
| 6B Biopsy controller interface (X204) | |

Fig. 10 Workstation main unit connectors of the Ultra 10


OPD00526

Ultra 1 Workstation

Connectors

CAUTION

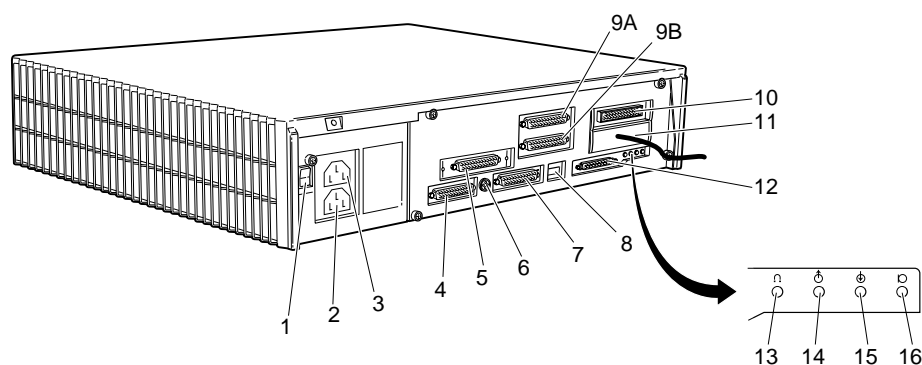Connecting the cable to the wrong connector may cause severe damage to the equipment.



- | | |
|--|--|
| 1 Power on / off | 10 Camera interface via biopsy controller (X203) |
| 2 Mains power supply | 11 Not used |
| 3 Do not use! | 12 SCSI interface for MO/CD-ROM drives |
| 4 Parallel port to local printer | 13 Not used |
| 5 Monitor interface | 14 Not used |
| 6 Keyboard interface | 15 Not used |
| 7 Not used | 16 Microphone |
| 8 Network connector | |
| 9A Mammomat interface via biopsy controller (X204) | |
| 9B Biopsy controller interface (X204) | |

Fig. 11 Workstation main unit connectors

Service and repair of the workstation is described in MAMMOMAT 3000 Opdima® Service Instructions Workstation (Ultra 1) SPB7-230.061.02.02.02.

OPD00624

CCD camera

Repair of the CCD camera is not allowed.

CAUTION

The CCD camera has to be handled with extreme care, it is very sensitive to mechanical shocks and temperature. When not connected, do not touch the pins in the camera contacts. In the camera shock and temperature sensors are integrated.

After installation of a new CCD camera the following has to be performed:

1. CCD camera calibration and maintenance, see CCD camera calibration and maintenance on Page 5 - 21.
2. Verifying the calibration of the biopsy unit, see Verifying the calibration of the biopsy unit on Page 10 - 1.
3. Check of Opdimas[®] AEC, see Check of Opdimas AEC on Page 10 - 2.
4. Check of resolution, see Check of resolution on Page 10 - 4.

Exchange of workstation

NOTICE

If a Ultra 1 is exchanged with a Ultra 10, notice that the trolley with part No. 64 30 921 for Ultra 1 does not fit Ultra 10. There is a new workstation table with part No. 65 34 106 for Ultra 10.

There is a kit available that consists of a Ultra 10 with workstation table, part No. 65 27 738. For installation of workstation table see Installation Instructions SPB7-230.031.18.01.02.

Procedure

1. Power off the workstation and MO unit. Disconnect all cables from the old workstation. Do not remove the monitor, keyboard or mouse.
2. Replace existing cable between the biopsy controller and the workstation with cable X204 part No. 64 98 138 (if the old workstation is a Ultra 1).
3. Place the Ultra 10 at a suitable place.
4. Connect cables to the Ultra 10.

NOTICE

Make sure that the SCSI cables are properly connected and fixed to the units and that the cable is not bent to much just behind the MO unit. Also, check that the MO unit is terminated correctly. If not, this might generate a lot of problems e.g. failure to write images to MO disk, corrupted MO disk or corrupted hard disk.

NOTICE

If the system is equipped with a combined CD/MO unit (Page 6 - 3), there is an extra MO shelf included in the Ultra 10/workstation table kit. The CD/MO unit shall be placed on that extra MO shelf and the MO shelf shall be mounted on the other side of the workstation table column. Adjust the height position of the shelf so it will not be squeezed in the table's lowest position (Fig. 12). If the CD/MO unit is installed to a Ultra 10, the internal CD unit of the Ultra 10 shall be used for future software upgrade.

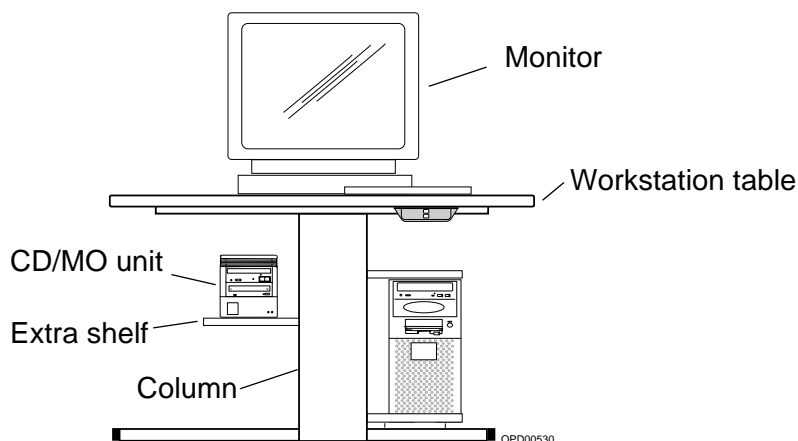


Fig. 12 Workstation table with an extra shelf

5. Switch on the MO unit and the workstation.
6. Log in as "poweroff" without password and wait until the OK prompt is shown on the screen.
7. Reboot the system by typing "`boot -r`" at the OK prompt.
8. Log in as "service".
9. Select Advanced service in the service dialog.
10. Select Software upgrade.
11. Restore hostname and data, see Page 6 - 24.

NOTICE

The local/backup MO disk used for restoring hostname must belong to the system being replaced if it shall be possible to use existing database including MO disks in the new workstation.

NOTICE

If the exchanged workstation had ASW 2.0/2.01 installed, steps 12.,13.,14. and 15. shall be performed if applicable.

If the exchanged workstation had ASW 2.1 installed, steps 13. and 15. shall be performed if applicable.

Steps 16. and 17. must be performed no matter what ASW the exchanged workstation had.

12. Restore Network setting, see Page 5 - 17.
13. Install DICOM option, see Page 5 - 16, a new license is needed due to new hardware. Apply and logout to make changes effective.
14. Restore DICOM nodes (see Supplement to the Instructions for Use MAMMOMAT 3000 - Opdimas[®], Page 78) by editing the IP addresses.
15. Restore printer configuration, see Page 5 - 24. Proceed with step 16.
16. Perform Unit tests, see Test of units on Page 5 - 12.
17. Perform Final procedures, see Final procedures on Page 6 - 24.

Reinstallation of software

After major failure, the software can be reinstalled. After this procedure no image data will be available on the temporary storage media. The system will prompt for appropriate local storage MO disks when viewing of images is requested.

For Ultra 1 only: If necessary, perform an installation of CD drive according to the Appendix 1.

Ultra 10 Pre installation

Installation of Operating Environment CD May 1999 or CD February 2000
(SUN CD distributed under Siemens part No. 65 27 696 (05/99) or 65 27 704 (02/00))

NOTICE

The system hostname (opdxxxx) is needed to perform the installation. Normally the hostname can be viewed on the screen when the system is switched on. If not, use the 4 digits from the serial No. labelled on the workstation main unit to create a dummy hostname. Restore the correct hostname from a MO disk from the original system after the installation is completed, see Restore hostname and data on Page 6 - 24.

NOTICE

Eject the MO disk and switch off the MO unit and the MAMMOMAT before performing the installation of the operating environment. This shall be done to avoid failure detection during installation of Solaris 2.5.1 HW 11/97.

1. Switch on the system.
2. Enter the keys:
Stop and A simultaneously.
This will bring you to the OK prompt.
3. Insert the Operating environment installation CD-ROM.
4. Type "boot cdrom" and press Enter.
5. In a few minutes, you will see a message asking you to select the language you want to see while running the installer.
Select English by selecting 1 and press Enter.
6. Would you like to use this installer?
Type "y" and press Enter.
7. Do you want to format/dev/dsk/c0t0d0?
Type "y" and press Enter.
8. Enter a swap partition size,
accept the default 512 MB and press Enter.

CAUTION

The next step will overwrite existing contents on the hard disk. Be certain that there is no data or software on the hard disk that you need to save.

9. Warning: All information on the disk will be erased:
Is this OK?
Type "y" and press Enter.

The system will now reboot and start the installation program, this will take approximately 5–10 minutes.

10. After the automatic reboot enter information about your system in response to questions in the System Identification windows. Step through the configuration of the operating system with the Next buttons and enter the data listed below.

- Host Name, set the original hostname (opdxxxx).

(If you do not know the hostname, use the 4 digits from the serial No. labelled on the workstation main unit to create a dummy hostname. Restore the correct hostname from a MO disk from the original system after the installation is completed, see Restore hostname and data on Page 6 - 24.)

-Network Connectivity, set No.

(The network is configured by the Opdima® Advanced service menus.)

-Time Zone, set the applicable geographic region.

- Date and Time, **Do not change** the default settings even though the settings are incorrect.

NOTICE

Changing the default settings for date and time may cause problems in the subsequent installation.

- Root Password, set no password by leaving the password fields empty. A warning message regarding password will appear, press OK.

- Power Management, accept the default values. (No power management.)

- Confirm Information by pressing confirm.

11. A welcome dialog will appear, press next.
12. A message will ask you which Solaris release you want to install. Select Solaris 2.5.1 Hardware: 11/97 and press next.
13. A summary dialog will appear, press next.

After you make this choice, the Operating Environment Installation CD will be ejected from the CD drive.

14. A message will ask you to insert Solaris 2.5.1 HW 11/97 CD (placed in a packing marked "Solaris 2.5.1 SOLMS-251WD999", included in the delivery. Insert the CD. Wait a few seconds and then press OK.

Installation of Solaris 2.5.1 HW 11/97

1. The message "Select Type of Install" will appear, select Custom Install and press next.
2. The message "System Locale Selection" will appear, select English (United States) (en_US) and press next.
3. The message "Select Solaris Cluster Configuration" will appear, select End User Solaris Software Group and press next.

4. The message "Lay Out File Systems" will appear, be sure that c0t0d0 (bootdisk) is applied in the column for Selected Disks and press next.
5. The message "Lay Out File Systems" will appear, select c0t0d0 and then press Modify.
6. The message "Disk c0t0d0" will appear. Remove the file system `/export/home`
Add all free memory to the `/` (root) file system.

Example with a 9 GB hard disk (figures can differ):

Disk c0t0d0 - 8692 MB

Slice	File system	Size
0	/	198
1	swap	512
2		
3		
4		
5		
6		
7	/export/home/	7982

Capacity: 8692 MB

Allocated: 8692 MB

Free: 0 MB

Rounding Error: 0 MB

Slice 7: Remove the file system `/export/home/` and size "7982".

Press Apply, 7982 is then added to "Free".

Slice 0: Edit the size of the `/` (root) file system from 198 to contain all memory available. In this example edit 198 to 8180 (198+7982). Observe that the value of the "Free" memory (7982) can differ.

Disk c0t0d0 - 8692 MB

Slice	File system	Size
0	/	8180
1	swap	512
2		
3		
4		
5		
6		
7		

Capacity: 8692 MB

Allocated: 8692 MB

Free: 0 MB

Rounding Error: 0 MB

Press Apply, press OK and press Next.

7. The message "Ready to Install" will appear, press Install Now.
The installation will take approximately 35–40 min.
Upon completion the CD-ROM will eject automatically.
Remove CD-ROM and close tray.
8. An installation summary will appear.
It is normal that the installation partially fails because the SUNWleoo package cannot be installed on Ultra 10. Select details to view the log file.
Press Details... for Solaris 2.5.1 HW 11/97.
Close to the end of the log file there is a summary of installed Solaris 2.5.1 packages. Verify that only SUNWleoo is listed under Solaris 2.5.1 packages partially installed. Press Exit to continue.
9. Press Reboot Now.

Installation of CDE

1. Login as "root".
No password is required.
2. Insert Desktop 1.1 CD-ROM from the media kit included in the delivery.
3. Type "`cd /cdrom/solaris_2_5_1_desktop_1_1/CDE/sparc`" and press Enter.
4. Type "`./install-cde`" and press Enter.
5. The message "Solaris CDE Installation Script Main Menu" will appear, select 1 (Default config) and press Enter.
6. The message "Begin CDE Installation now?" will appear, select Yes and press Enter.
This will take approximately 3–5 min.
7. To eject the Desktop 1.1 CD-ROM, type "`cd /`" and press Enter. Type "`eject`" and press Enter.
Remove the CD-ROM from the CD drive.
8. Switch on the MO unit before rebooting the system. Do not insert the MO disk.
9. Type "`reboot`" and press Enter.
10. Proceed the installation with **Ultra 10/Ultra 1 Installation of Opdima ASW 2.1 (part No. 65 27 639)** on Page 6 - 22.

Ultra 1 Pre installation

Installation of Opdima ASW 1.3

NOTICE

The system hostname (opdxxxx) is needed to perform the installation. Normally the hostname can be viewed on the screen when the system is switched on. If not, use the 4 digits from the serial No. labelled on the workstation main unit to create a dummy hostname. Restore the correct hostname from a MO disk from the original system after the installation is completed, see Restore hostname and data on Page 6 - 24.

NOTICE

Eject the MO disk to avoid loss of information.

1. Connect and install the CD drive if the system has an separate external one, see Appendix 1.
2. Switch on the system.
3. Enter the keys:
Stop and A simultaneously.
This will bring you to the OK prompt.
4. Insert the ASW 1.3 CD-ROM.
5. Type "boot cdrom" and press Enter.
The initial boot phase takes about 3 to 4 minutes.

CAUTION

The next step will overwrite existing contents on the hard disk. Be certain that there is no data or software on the hard disk that you need to save.

6. Confirm formatting of the hard disk by typing "y"(es) and press Enter.
7. If ASW 2.1 has been previously installed on the system the hostname is automatically retrieved from the NVRAM. Otherwise, the Opdima® system may tell you to Please enter the 4-digit site number. Enter the 4 digits from the hostname and confirm the hostname: opdxxxx (where xxxx=site number) by pressing "y"(es) and Enter.
(If you do not know the hostname, use the 4 digits from the serial No. labelled on the workstation main unit to create a dummy hostname. Restore the correct hostname from a MO disk from the original system after the installation is completed, see Restore hostname and data on Page 6 - 24.)
8. The installation process will start and take about 10 minutes.
9. Specify system type as stand alone by typing "s" and press Enter.

NOTICE

Network settings is performed in Advanced service.

10. Confirm the configuration by pressing “y”(es) and press Enter.
The system will now reboot and the login dialog will automatically be displayed when the configuration is ready.

NOTICE

Remove the CD-ROM during reboot.

11. Proceed the installation with **Ultra 10/Ultra 1 Installation of Opdima ASW 2.1 (part No. 65 27 639)** on Page 6 - 22.

Ultra 10/Ultra 1 Installation of Opdima ASW 2.1 (part No. 65 27 639)

NOTICE

If there is an existing patient database on the system, it is very important to do Restore hostname and data on Page 6 - 24 directly after the installation of ASW 2.1 is completed. If not, the existing database will be lost.

1. Login as "root".
For Ultra 10 only: No password is required.
For Ultra 1 only: The root password is "v1armkorv".

NOTICE

Make sure that the Num Lock and Caps Lock are deactivated.

2. Place the mouse pointer on the screen background and press the rightmost mouse button.
3. Select Programs and Terminal.
A window with a UNIX shell prompt will appear.
4. Insert the Opdima® ASW 2.1 CD, and wait a few seconds.
5. Change directory by typing
"cd /cdrom/opdima.2.1_6527639" and press Enter.
6. Start the installation script by typing ". /install" and press Enter.
7. The message "Do you wish to install the required patches?" will appear, answer with yes by typing "y" and press Enter.
For Ultra 10 only: This will take approximately 40 minutes.
For Ultra 1 only: This will take approximately 5 minutes.
8. The message "To continue the installation, you will have to reboot. Reboot now?" will appear, answer with yes by typing "y" and press Enter.
9. Press Enter to continue the installation.
10. After reboot, repeat steps 1. to 3. and 5. to 6.
11. The message "Do you wish to install the Opdima® software and misc. binaries?" will appear, answer with yes by typing "y" and press Enter.
Loading software from CD-ROM will take approximately 5 minutes.

NOTICE

Check that the packages "OPDcore", "OPDlang", "OPDmisc", "OPDmodisc", "OPDsol251", and for Ultra 10 "OPDu-5-10" and for Ultra 1 "OPDu-1" are completely installed, which will be displayed on the system monitor.

12. The message "To complete the installation, you will have to reboot. Reboot now?" will appear, answer with yes by typing "y" and press Enter.
13. Remove the ASW 2.1 CD-ROM and press Enter to continue.

⚠ WARNING

Do not enter stereo/spot or database mode before Restore Disk has been performed to avoid loss of database information.
Be sure that no MO disk is inserted in the MO unit.

14. The system reboots.
If the system has an external CD drive, power off the system and disconnect the CD drive after the software reinstallation.
15. Proceed with Restore hostname and data on Page 6 - 24.
16. Proceed with Final procedures on Page 6 - 24.

Restore hostname and data

1. Log in as service user.
2. Select Service and Advanced service.
3. Insert the most recently used local/backup MO disk in the MO unit.
4. If the 4 digits from the serial No. labelled on the workstation main unit has been used as hostname and there are doubts if that is the correct hostname (which is the case when the computer has been exchanged), select Software upgrade and Restore hostname.

NOTICE

Restore hostname from most recently used local/backup MO disk.

5. Select Restore Disk and press Enter to restore data from MO disk (see Page 5 - 20). This copies the database from the MO disk to the hard disk. It is possible to use a local storage MO disk or a backup MO disk.

WARNING

When restoring data, be sure to insert the MO disk that was used most recently before the reinstallation of software.

Date for last update of MO disk is displayed on the monitor before confirmation of restore.

If the system is restored from an older MO disk, the most recent folders will be lost from the database and the numbering of the new MO disks can be incorrect. Please contact HSC for more information.

6. Select Restore and press Enter to copy the data to the hard disk.

NOTICE

Set-up the printer once again after software reinstallation. Printer information can not be restored.

Final procedures

1. Check in Database mode that the latest examination performed is available and the images can be viewed.
2. Check camera image quality in normal and high resolution. Use calibration exposure data, see Appendix 3.
3. If the image shows defects e.g. lines, spots or irregularities such as 45 degree-mesh pattern, calibrate the camera according to CCD camera calibration and maintenance on Page 5 - 21.
4. Check biopsy accuracy with the test phantom, if the accuracy is not within $\pm 1\text{mm}$ in (x, y, z), calibrate the biopsy unit according to Calibration of the biopsy unit on Page 5 - 3.

Workstation

The system displays error and warning messages on the screen for abnormal situations:

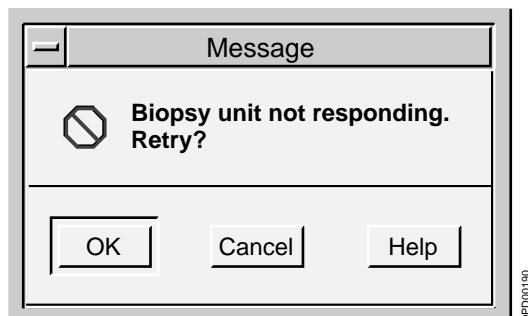


Fig. 1 Error message dialog

The table below shows some of the messages that can occur on the workstation screen. For other messages, see Messages in Supplement to the Instructions for Use MAMMOMAT 3000 - Opdima®.

Message	Possible cause and action to take
An error occurred while calibrating. Please press EXIT and try again!	Failure during stereo calibration. Press Exit and try again.
Backup could not be performed. Press OK to continue.	Failure during the backup procedure. Press OK. It could be problems with cabling, SCSI termination or dust in the MO unit. Check power and cabling and retry.
Biopsy unit not responding. Retry?	No contact with biopsy unit. Retry.
Camera exposure timeout.	The time between pressing Acquire from the workstation and releasing the exposure from the MAMMOMAT generator has exceeded 60 seconds.
Camera not responding. Retry?	No contact with camera. Retry.
Cannot use biopsy unit.	Communication error or hardware failure in biopsy unit. Check power and cabling and retry.
Cannot use camera.	Communication error or hardware failure in CCD camera. Check power and cabling and retry.
Cannot use Mammomat.	Hardware failure in MAMMOMAT. Check power and cabling and retry.
Double exposure with improper stereo angles!	A stereo pair has, possibly by mistake, been exposed with the same stereo angle. If decided, acquire and expose again with the correct stereo angles.
Error transmitting targets to the biopsy unit.	The system has detected an error while transmitting target to the biopsy unit. Check connections. Try again.

Message	Possible cause and action to take
Exposure aborted.	The exposure failed due to MAMMOMAT or AEC error.
Exposure values out of range.	The values given by the MAMMOMAT are out of range. Possible failure in MAMMOMAT or installation.
Failed to calibrate biopsy unit.	Failure during calibration. Press Exit and try again.
Failed to copy folder(s).	The target disk may be full or corrupt.
Failed to initialize for Network Setup.	Failure during Network Setup. Press OK. Check power and cabling and retry.
Failed to initiate biopsy unit.	Communication error at start-up. Check connections to the biopsy controller and retry.
Failed to store on MO disk.	Check cabling, MO unit and MO disk. MO disk may be full.
Failed to transmit target.	The system has detected an error while transmitting target. Check connections to the biopsy controller and retry.
Find and adjust the reference marks, then press OK.	Exposure has been performed, the reference marks cannot be found. Adjust contrast and brightness to visualize them. Click with the pointer on the reference T and press OK.
Internal error.	A software failure has occurred. Restart the system.
Internal error. Retry?	A software failure has occurred. Retry.
Mammomat not responding. Retry?	No contact during start-up. Retry.
Storage MO disk not available.	Check cabling, MO unit and MO disk. Check if MO disk is inserted.
System Error.	A system failure has occurred. Restart the system.
Warning: Camera installation not complete	Incorrect camera calibration. Perform camera calibration and maintenance, see Page 5 - 21.
XX images expected, YY found.	The system has detected a failure in the database. If backup exists, restore backup data, see Restoring backup data in Supplement to the Instructions for Use MAMMOMAT 3000 - Opdima®

Control panel

The MAMMOMAT generator control panel gives the user additional information at abnormal situations.

Error message	Possible cause and action to take
Error message 777: "Exposure sequence is aborted by Opdima®".	<ul style="list-style-type: none"> - The time between pressing Acquire from the workstation and releasing the exposure from the MAMMOMAT generator has exceeded 60 seconds. Press the Limit button on MAMMOMAT generator control panel to continue. - The exposure parameters (kV, mAs) chosen will give a too weak beam quality for the object. Press the Limit button on MAMMOMAT generator control panel to continue. - The exposure parameters (kV, mAs) chosen will give a too hard beam quality for the object. E.g. if AEC mode is used and the estimated mAs is lower than 5 mAs. May happen when using the stereo calibration phantom. Press the Limit button on MAMMOMAT generator control panel to continue and choose lowest kV in AEC mode or set manual exposure parameters. - The Cancel button has been pressed while "Expose at MAMMOMAT" or "Exposure preparation started" were displayed. Press the Limit button on MAMMOMAT generator control panel to continue.
Error message 778: "No communication present between MAMMOMAT and Opdima®".	Press Retry from workstation or press the Limit button on MAMMOMAT generator control panel to continue.
Error message 779: "The mAs value for the main exposure in Opdima® AEC has not been received in time".	Communication problem. Press the Limit button on MAMMOMAT generator control panel to continue.
External diaphragm lamp is lit on the generator control panel.	Opdima® external diaphragm not mounted. Mount Opdima® external diaphragm.

WARNING

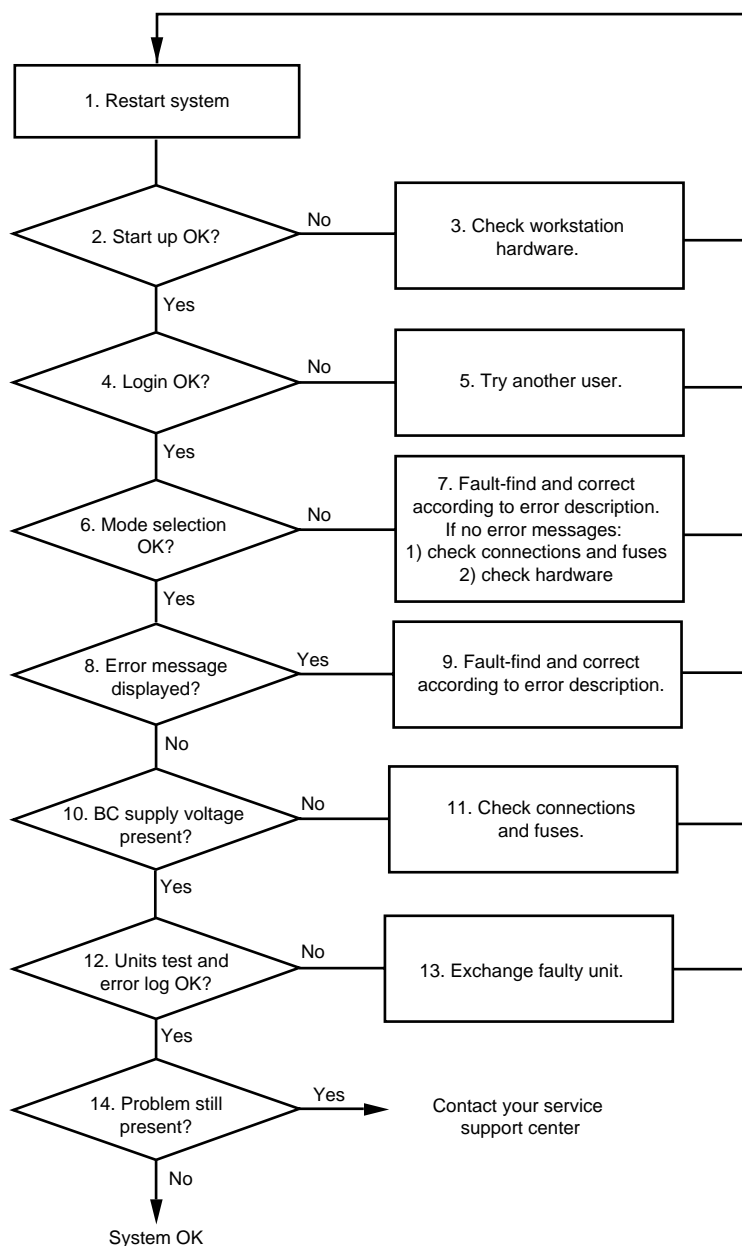
When resetting error 778 with the limit button on the MAMMOMAT, the MAMMOMAT will switch to film exposure mode ("OP" disappears from the display on the generator). To resume Opdima® mode; in the control dialog of the Opdima® system, select Exit and choose the required mode (Stereo or Spot) in the Mode selection dialog.

This page intentionally left blank.

Fault isolation chart

NOTICE

This chart is to be regarded as a guide only, and not as a complete service tool.



OPD00525

Fig. 1 Fault isolation chart

Explanations

Step	Explanation
1	Restart the system.
2	The workstation hardware is tested automatically when booting the workstation, (the Opdima® hardware is not tested). If an error is detected, the type of error is displayed on the screen.
3	Check the connectors of the hardware units. Exchange faulty units.
4	Log in as service.
5	The user name and password used could have been invalid.
6	Depending on selected mode, tests are performed automatically.
7	Take measures according to the displayed error message. If no error messages: check connections and fuses, check hardware.
8	Error messages are displayed on the workstation screen if the system detects an error.
9	Take measures according to the displayed error message.
10	Check LEDs indicating voltage present in the biopsy controller.
11	Check connections and fuses in the biopsy controller.
12	Perform units test from the workstation service mode. Check the log for abnormal occurrences.
13	Change faulty unit according to instructions.
14	If the problem is still present, contact your service support center.

General

In general, always check the log files for more information, see Log inspection on Page 5 - 28.

NOTICE

The MO disks used in the Opdima® system are equipped with a mechanical write protection lock. This lock shall never be used since the Opdima® system always attempts to mount the MO disk as a read/write device.

If the write protection lock is activated, the behavior is undefined. There may be a lot of error messages written across the screen, hiding the Opdima® dialogs (windows), internal error may occur and/or the system may reboot spontaneously.

If error messages occur on top of the normal Opdima® screen, wait until no more error messages are written and the activity lamp on the MO unit is turned off. This may take a while (30 minutes) depending on the operation that initiated the error messages. Then eject the MO disk if possible, otherwise see below.

All effort shall be taken to try to find the control dialog (window) and Exit the normal way. The dialogs may be hidden by error messages, but can be found by pointing and clicking on screen. Log out and log in as "poweroff" to shut down the system. If necessary, power off the MO unit and turn it back on to be able to eject the disk. Then turn on the workstation again.

Biopsy unit

Biopsy unit not responding

- Check power switch on biopsy controller.
- No contact during start-up. Retry.
- Check cables.
- Check biopsy controller, fuses and power.
- Replace board D200. You will always get Biopsy unit not responding when replacing D200. LedV2 NVM error will light up, goes out after calibration.
- Recalibrate biopsy unit.
- If “Biopsy unit not responding” occurs after calibration, proceed with retry or log in as “user”.

Problem with calibration

- Check that you have all targets in the image, use 0.7 magnification.
- Check that you use the right needle length.
- Check that you have set the same target both on the biopsy unit and in the workstation.

Camera

Camera not responding

- No contact during start-up. Retry.
- Check that MAMMOMAT and Biopsy controller are powered on.
- Power off and reboot Opdima® workstation.
- **For Ultra 10 only:** View error log, if message cannot reach CCD camera on/dev/CCD0 appears, check that D101 is properly connected inside the workstation.
- Check cables.
- Check biopsy controller, fuses and power.
- Change D100/D101.
- Change D200.

Cannot use camera

- Hardware failure in CCD camera. Retry.
- Check that MAMMOMAT and Biopsy controller are powered on.
- Power off and reboot Opdima® workstation.
- Check cables.
- Check that board D100/D101 is properly connected.
- Change D100/D101.
- Change camera.

Workstation

If “Bogus file system” appears

- Check cables, see SCSI cables, CD-ROM and Cables general in MAMMOMAT 1000/3000/3000 Nova - Opdima® Maintenance Instructions.
- Try to reboot.
- Reload the software,
for Ultra 10, see Page 6 - 16 and for Ultra 1, see Page 6 - 20.

Problem with date in database mode

- Check the Country Settings Date Pattern: e.g. m/d/y. This defines only the order that the date will be in. Do not type mm/dd/yy.

No images displayed on monitor

- If no image on the monitor, check that input select is set to 1.

For Ultra 10 only:

NOTICE

If removing the cover of the workstation main unit, check that no copper grounding clip comes loose.

Tests

- Power on failure. Ultra 10 Service Manual, chapter 4.
- Video output failure. Ultra 10 Service Manual, chapter 4.
- Hard drive or CD drive failure. Ultra 10 Service Manual, chapter 4.
Comment: When replacing the hard drive, reload the complete software.
- Power supply test. Ultra 10 Service Manual, chapter 4.

For more information in general about Ultra 10 see Ultra 10 Service Manual.

Ultra 1 If the workstation does not start after an unplanned power cut off

- Unplug the AC power cord from the wall socket.
- Let it be unplugged for 20 minutes.
- Connect the AC power cord and switch on the workstation.

Image Quality

Problems during calibration

- Check that gridlines are not visible in images taken with a table without grid e.g. biopsy unit. If visible, see CCD camera calibration and maintenance on Page 5 - 21.
- Calibration not completed takes unexpectedly long time. This can be a defective camera or a mechanical alignment problem. Try to calibrate without external collimator, activate switch or use a collimator with larger hole e.g. Axilla collimator.

NOTICE

Always use Opdima® collimator in normal use.

Quadrant difference on patient images

When doing images on homogeneous acrylic plastic (PMMA) plates there can be a slight difference between the quadrants, which will not be visible in patient images. Make sure that quadrant difference exists in patient images before changing the CCD camera.

- Calibrate the CCD camera.
- If it is not disappearing, please contact the headquarter support center, HSC 24/USC, for investigation.

Quadrant missing

- Try to calibrate the CCD camera.
- Please contact the headquarter support center, HSC 24/USC, for investigation.

White line/dot

- Calibrate the CCD camera, move acrylic plastic plates after each calibration exposure.
- If it is not disappearing, please contact the headquarter support center, HSC 24/USC, for investigation.

MAMMOMAT**Cannot use MAMMOMAT**

- Try again.
- Switch off/on the MAMMOMAT.
- Power off and reboot the Opdima® workstation.
- Check biopsy controller, fuses, and power.
- Check cables.
- Check D707, must be installed.

Network problem

DICOM (option), checklist:

- A RJ-45 network connector is needed.
- If the distance between the Opdima® workstation and the network connector is longer than 3,5 m (3,5 m cable included in the package) a UTP-5 cable of maximum 50 m is needed.
- The license key to enable the DICOM (option) is written on the invoice and the dispatch note (if not, please contact Siemens-Eléma AB mammography department in Solna, Sweden for further information).
- To be able to receive the optional DICOM license key the mammography logistics department will request for the serial No. and the host ID for the workstation in question. How to receive these parameters is described in Instruction for obtaining the systems host ID and hostname on Page 5 - 15.
- The Opdima® system needs the following information to be connected to the network:
IP address: (e.g. 123.123.123.11)
Netmask: (e.g. 255.255.255.0)
Default router: (e.g. 123.123.123.22)

NOTICE

The IP address and the default router may not contain any initial 0, e.g. 10.10.10.10 is OK but 010.010.010.01 is not OK.

For further information see DICOM Conformance Statement for Opdima® ASW 2.0 which is available on Intranet.

- Network information about the DICOM nodes, which will be connected to the Opdima® system.
The following information for all DICOM nodes is necessary:
Hostname: (e.g. opdxxxx)
IP address: (e.g. 123.123.123.11)
DICOM Port: (e.g. 1984 for Opdima)
DICOM AET: (e.g. opdxxxx_SRV)

NOTICE

For Opdima®, hostname, DICOM Port and DICOM AET are fixed.

xxxx is the serial No. for the workstation. The serial No. is labelled on the cover of the main unit (last four digits). The last four digits of the hostname normally consists of the last four digits of the serial No. This will not be the case if there has been a hostname change due to a workstation exchange.

General Information about networking

- Project Manager Handbook: TD home page, <http://www-td.med.siemens.de>
Product Information
Planning
Project Manager Handbook

Printer

Printer connected directly to Opdima

Print problem:

- Check the crossover network cable.
- Check IP address

Printer connected to network

Print problem:

- Check cables, normal network cable.
- IP address from network administrator at the hospital.
- When defining the printer in advanced service use PS for all postscript printers and XWD for Codonics printer.

MO unit / CD drive

- Check that correct media is used, see spare part list.
- Check that write protection is not activated.

All problems with MO disk

- Identify MO disk, see Page 5 - 19.
If more than one disk type or disk number is reported please contact the headquarter support center, HSC 24/USC, for investigation.
- Check and, if necessary, repair MO disk, see Page 5 - 19.
- Make sure that the SCSI cables are properly connected and fixed to the units and that the cable is not bent to much just behind the MO unit.
- Check that the MO unit is terminated correctly.

If the message “Failed to store on MO disk” appears

- Make sure that local storage disk is inserted.
- Make sure that MO unit is switched on.
- Identify MO disk, see Page 5 - 19.
- Check if MO disk is full.
- Check and, if necessary, repair MO disk, see Page 5 - 19.
- Try to save on a new MO disk.
- Power off and reboot the Opdima® workstation.
- Check cables, termination. Change MO unit.

If the message “Cannot read disk Prepare disk?” appears

- If it is a new disk, press OK (this will format the disk). Else, press Cancel.
- Identify MO disk, see Page 5 - 19.
- Check and, if necessary, repair MO disk, see Page 5 - 19.
- Power off and reboot the Opdima® workstation.
- Check cables and termination, see MAMMOMAT 1000/3000/3000 Nova - Opdima® Maintenance Instructions.
- Change MO unit.

Cannot communicate with MO unit

Test of SCSI bus



WARNING

Test of SCSI bus can destroy the hard disk content. The software may have to be reinstalled.

- Press the buttons Stop and A simultaneously as soon as possible when the system is booting up, you will then get the OK prompt.



If Stop A is not pressed directly when the system starts to boot, the probe-scsi-all command will hang the system without showing any SCSI-devices. Power supply switch must be used to turn off the system.

- Type “probe-scsi-all”, the system should show all the drives connected to the bus.

For Ultra 10 only:

MO unit	Target 1
(The hard disk and CD drive are connected to an internal ide bus)	

For Ultra 1 only:

Hard disk	Target 0
Mo unit	Target 1
CD drive	Target 6

Verifying the calibration of the biopsy unit

Perform stereo examination with stereo calibration phantom.

1. Prepare the system according to Setting up the system for stereo examination in Supplement to the Instructions for Use MAMMOMAT 3000 - Opdima®.
2. Perform step 1. to 25. of a stereo examination on the stereo calibration phantom (Fig. 4, Page 5 - 4), see Performing stereo examination in Supplement to the Instructions for Use MAMMOMAT 3000 -Opdima®.

Use a fictive patient ID. Record the patient ID in Appendix 3.

3. Adjust the adjustment knobs on the biopsy unit so that the needle will not hit the stereo calibration phantom. Insert fine needle.
4. Use the adjustment knobs on the biopsy unit to adjust the needle to the stereo calibration phantom target. The needle shall be positioned according to Fig. 5, Page 5 - 4.
5. Check that the biopsy unit displays do not show values more than ± 1 mm. If they do, perform a biopsy calibration. After calibration of the biopsy unit, perform step 1. to 5. again. Fill in your values in the test protocol for repetitiveness (see Appendix 3).

Check of Opdima AEC

The objective is to assess the performance of the AEC (automatic exposure control) function with different exposure conditions.

In this test, images shall be obtained using different combinations of kVp, anode/filter material and object thickness. The AEC is designed to control the exposure so that the mean AD value of the images obtained of homogeneous phantoms shall be close to 800 independently of these parameters. For some combinations of kVp, object thickness and anode/filter material, the minimum mAs value (12 mAs) or maximum mAs value (560 mAs, or lower due to MAMMOMAT limits) is the limiting factor for proper AEC functionality. These combinations are of no clinical relevance and should not be tested, and have been marked with a gray color in the test report tables (Appendix 3).

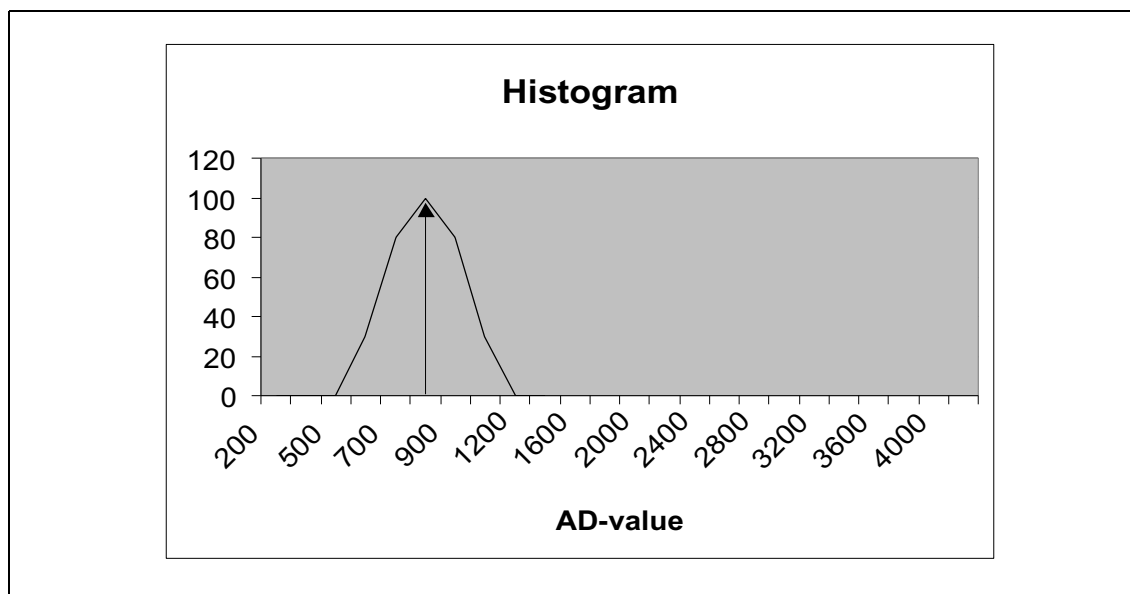


Fig. 1 Histogram

Procedure

1. Start a new examination.
2. Enter Spot mode.
3. Use the table that will be used by the customer, images shall be taken at normal resolution, and the mammography unit shall be used with a grid object table.
4. Initially, place 20 mm of PMMA on the object table.
5. The first image shall be obtained using 26 kV and Mo/Mo as exposure parameters.
6. Examine the image and measure the mean value. The gray scale parameters (center and width) are shown in the upper bar of the Opdima® window. The center value is approximately the mean value and shall be used as such in this test. Note that these values (gray scale center and width) are automatically set to optimal when an image is captured, but will be affected by adjustments made to brightness or contrast in Opdima®. Because of this, no adjustments shall be done to these parameters before reading the AD center value, in order to get a correct reading.

7. Fill in the mean (center) value in the appropriate space in the test report tables (Appendix 3).
8. Repeat step 4. to 7. with all combinations of object thickness, kV parameter and anode/filter combination that are found in the test report tables in the Appendix 3.

Performance Criteria

The measured mean value shall be equal to $800 \pm 15\%$ (above 680 and below 920) for all examined combinations of object table, anode/filter material, kV value and object thickness.

NOTICE

If the exposure times are long it can be possible to distinguish that the image is divided into a pattern of “blocks” see Fig. 2. The system is calibrated for clinical exposure times which shall be kept short in order to avoid movement blurredness of the image and also to avoid high dose exposure of the patient.

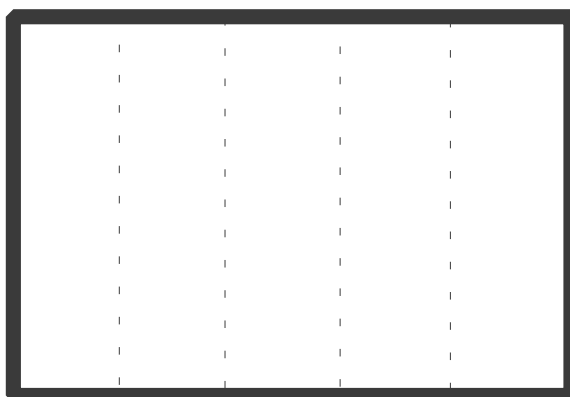


Fig. 2 Pattern of “blocks”

Check of resolution

The objective is to ensure that a sufficient spatial resolution is obtained with the Opdima® system.

Procedure

1. The grid object table shall be used.
2. Select Mo/Mo.
3. To get a proper mAs-value, perform a test exposure in AEC-mode with 4.5 cm PMMA without the bar pattern (it might interfere with the sensitive area of the AEC). Use a kVp that is commonly used clinically and record this value.
4. Manually select (on the control panel of the MAMMOMAT) the mAs-value that is closest to the measured AEC mAs value and use it during the following measurements (with AEC mode disabled).
5. Position the “Bar-pattern” on 4.5 cm of PMMA, 1 cm from the chest wall edge. Position the bar pattern diagonally, i.e. with the bars at an angle approximately 45 degrees to the anode-cathode axis.¹ Make an exposure at normal resolution.
6. Use the magnification tool and change grayscale settings to optimize viewing conditions for this evaluation. The grayscale settings are changed using the “brightness” and “contrast” controls in Opdima®, either using the scroll bars at the right side of the screen, or by clicking with the middle button of the mouse and dragging. Adjust these setting in order to get the best possible image of the bar pattern. For example, if using 27 kVp and Normal resolution, a gray scale center of 380 and of width 110 might be useful. Note that these values are only guidelines - depending on the specific systems, other values might be more appropriate.
7. Examine the image and determine the highest line-pair resolution where dark/light lines can still be differentiated, and note this value in the test report tables (Appendix 3).
8. Repeat steps 3. to 7. with a W/Rh² target/filter combination, if used clinically.
9. Repeat steps 2. to 6., but make the exposure at high resolution this time.

Performance Criteria

The measured resolution shall be at least 10 lp/mm in normal mode, and at least 13 lp/mm in high-resolution mode³.

1. This test is designed to ensure the performance of the Opdima system and not the MAMMOMAT, which should be tested separately. Therefore, the test is performed with the bar-pattern positioned at a single angle. A test of the bar-pattern positioned diagonally is sufficient to ensure quality of the Opdima components, and differences in resolution at different angles (parallel vs. perpendicular to anode-cathode axis) is due to MAMMOMAT limitations, and shall therefore be tested and verified in the QC of the MAMMOMAT. Interference between the bar-pattern and the detector grids might affect measurements if they are performed with parallel or perpendicular orientation.

2. It is not necessary to perform the test with a M/Rh-combination since different filters do not have an impact on resolution.

3. The maximum theoretical resolution for the detector in Normal Resolution mode is 12.4 lp/mm, and the detector (part of Opdima) will be the limiting factor. For High Resolution mode the maximum theoretical resolution is > 20 lp/mm, and the x-ray tube (part of the MAMMOMAT) will be the limiting factor.

Protective earth measurement

Protective earth measurement shall be performed after completed hardware service. Record type, serial number and calibration date of the instruments together with measured values in the service report.

Biopsy unit

Use the protective ground wire tester ($\leq 6V/10-25A$).

1. Measure the resistance between the ground connection of the stand and the measuring points on the biopsy unit according to Fig. 3. Use a sharp-pointed tool to penetrate the finishing paint. Fill in your values in the test protocol for protective earth measurement (see Appendix 3).

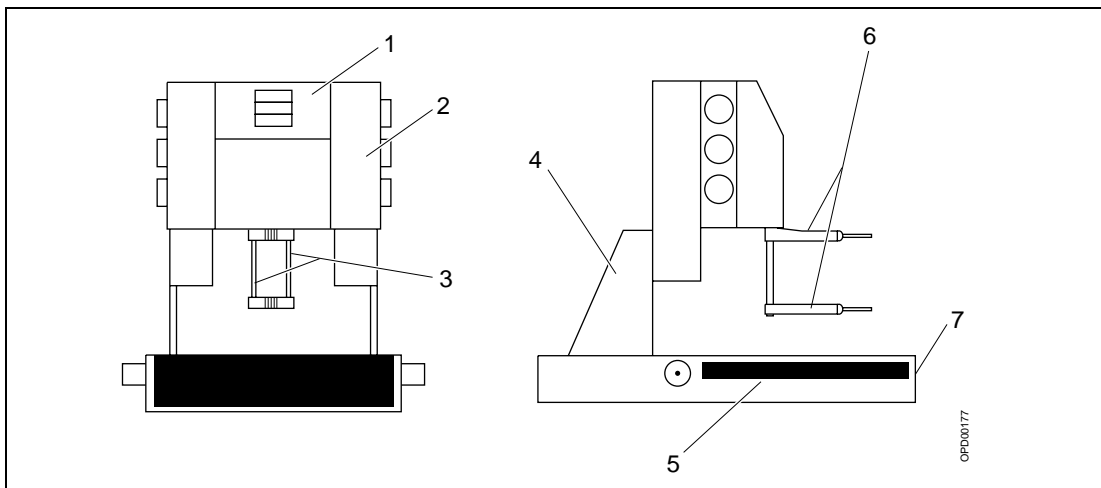


Fig. 3 Measuring points on the biopsy unit

The value of the resistance measured must not exceed 0.1Ω , except for measuring points 3 and 6 for which no limit is required (the test performed complies with IEC 601-1 subclause 18A).

Biopsy controller

Use the protective ground wire tester ($\leq 6V/10-25A$).

1. Measure the resistance between the ground connection of the stand and the cover of the biopsy controller. Use a sharp-pointed tool to penetrate the finishing paint. Fill in your values in the test protocol for protective earth measurement (see Appendix 3).

The value of the resistance measured must not exceed 0.1Ω .

CCD camera

Use the protective ground wire tester ($\leq 6V/10-25A$).

1. Measure the resistance between the ground connection of the stand and measuring points on the CCD camera according to Fig. 4. Use a sharp-pointed

tool to penetrate the finishing paint. Fill in your values in the test protocol for protective earth measurement (see Appendix 3).

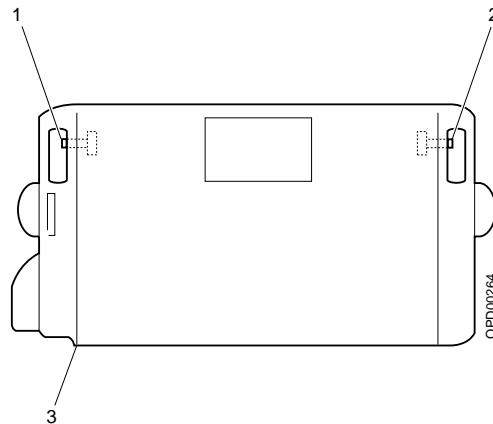


Fig. 4 Measuring points on the biopsy unit

The value of the resistance measured must not exceed 0.1 Ω .

Biopsy controller cable duct

Use the protective ground wire tester ($\leq 6\text{V}/10\text{--}25\text{A}$).

1. Measure the resistance between the ground connection of the stand and the cover of the biopsy controller cable duct. Use a sharp-pointed tool to penetrate the finishing paint. Fill in your values in the test protocol for protective earth measurement (see Appendix 3).

The value of the resistance measured must not exceed 0.1 Ω .

Chap.	Name of chapter	U1	U10	Both	Changes to previous version
All					Service Instructions for Ultra 1 and Ultra 10 have been put together in this document.
0				X	Completely revised.
1	General			X	The document is valid for ASW 2.1.
	Training of customer support engineer			X	Training of customer support engineer is new.
	Documents required			X	Maintenance Instructions is a required document.
	Documents required	X			Installation and Start-Up Instructions is not a required document.
	Ultra 10 CD-ROMs required		X		<ul style="list-style-type: none"> - There are two new operating environment installation CDs. - CD-ROM ASW 2.01 is replaced with CD-ROM ASW 2.1. - CD-ROM serial port configuration patch is removed.
	Ultra 1 CD-ROMs required	X			<ul style="list-style-type: none"> - CD-ROM ASW 2.0/2.01 is replaced with CD-ROM ASW 2.1. - CD-ROM serial port configuration patch is removed.
2	Workstation	X			CD drive is included in the workstation.
	Workstation			X	Text change.
5	Calibration of the biopsy unit			X	<ul style="list-style-type: none"> - Two new notices. - Figure 4 has been changed.
	Backup of temporary storage media		X		Text change.
	Country settings			X	<ul style="list-style-type: none"> - There are four languages available in Country settings. - Time settings is new in Country settings. - Text change in table.
	Country settings	X			Institution name and Department name are new.
	Test of units	X			The notice is new.
	Software upgrade			X	<ul style="list-style-type: none"> - Restore hostname is new in Software upgrade. - Update eeprom is new in Software upgrade. - Add language is removed from Software upgrade. - Text change in table.
	New software version			X	<ul style="list-style-type: none"> - Figure has been changed. - Text change in caution.

	Instruction for obtaining the systems host ID and hostname			X	Instruction for obtaining the systems host ID and hostname has been moved here from Installation and Start-Up instructions and is new for Ultra 1.
	DICOM option			X	- Dicom option has been moved here from Installation and Start-Up Instructions. - Text change.
	Network settings			X	- Network settings has been moved here from Installation and Start-Up Instructions. - Text change.
	Miscellaneous			X	New options in Miscellaneous: Inspect database contents and MO disk utilities.
	Inspect database contents			X	Inspect database contents is new
	MO disk utilities			X	MO disk utilities is new.
	Restoring data from MO disk			X	Text change in warning.
	CCD camera calibration and maintenance			X	Text change.
	Disk cache settings			X	Text change in table.
	Use of printer setup			X	- Use of printer setup has been moved here from Installation and Start-Up Instructions. - Text change.
	Log administration			X	Text change.
6	General			X	- Caution is new. - Text change.
	MO unit (no CD)	X			MO unit is new.
	Removing the MO unit			X	- Caution is removed. - Notice about C-MOS guidelines is removed. - Text change.
	Installing the MO unit			X	- Notice is new. - Figure has been changed. - Text change in table and in text.
	CD/MO unit	X			CD/MO unit is new.
	Removing the CD/MO unit	X			- Caution and first notice are removed. - Text change in notice.
	Installing the CD/MO unit	X			- First notice is new. - Text change in second notice and in table. - Figure has been changed.
	Replacement of board D200			X	Notice is new.
	Ultra 10 workstation		X		Second caution is new.

	Ultra 10 workstation			X	Figure has been changed.
	Ultra 1 workstation	X			Figure has been changed.
	Exchange of workstation			X	Exchange of workstation is new.
	Reinstallation of software			X	Text change.
	Restore hostname and data			X	Restore hostname and data is new.
	Final procedures			X	Text changes in Final procedures.
7	Messages			X	- Messages have been moved here from chapter 8. - Text changes in table.
8	Fault isolation chart			X	- Fault isolation chart has been moved here from chapter 5. - Text change in tables.
9	Troubleshooting guide			X	Troubleshooting guide is new.
10	Verifying the calibration of the biopsy unit			X	Text change.
	Check of Opdim AEC			X	Text change.
	Check of resolution.			X	Text change.
A1	Appendix 1	X			Appendix 1 is new.
A2	Appendix 2			X	Appendix 2 is new.
A3	Appendix 3			X	- Protocol for CCD camera calibration is new. - Protocols for image quality have been changed.

This page intentionally left blank.

Ultra 1 Installation of CD drive, part no 65 38 180

NOTICE

The CD drive should only be installed during software installation.

NOTICE

When the CD drive is in use, move the external termination from the MO unit to the CD drive.

Dip switch settings for CD drive, part no 65 38 180:

No. 1	ON	Parity check
No. 2	OFF	Internal termination shall not be used
No. 3	OFF	Test (factory use only)
No. 4	ON	Block size (must be ON in an UNIX environment)
No. 5.	OFF	Eject lock

1. System off.
2. Connect the SCSI cable from the CD drive according to Fig. 2.
3. Move the external termination from the Mo unit to the CD drive.
4. Set SCSI ID 6 for the CD drive and set the dip switches according to the table above.
5. Switch on the CD drive and the MO unit.
6. Switch on the workstation by pressing the power on key.

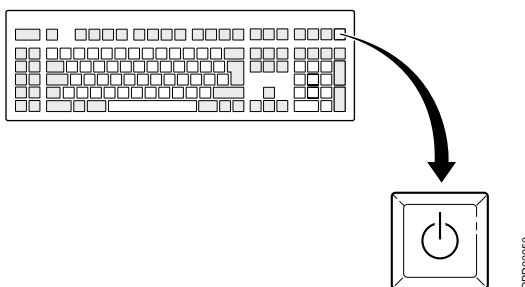


Fig. 1 Power on key

NOTICE

If the system does not find the CD drive, reboot the system by typing "boot -r" at the OK prompt.

NOTICE

Put the external termination back to the MO unit when the CD drive is disconnected.

Ultra 1 Cables

Use the following SCSI cables between the CD drive and the MO unit/workstation.

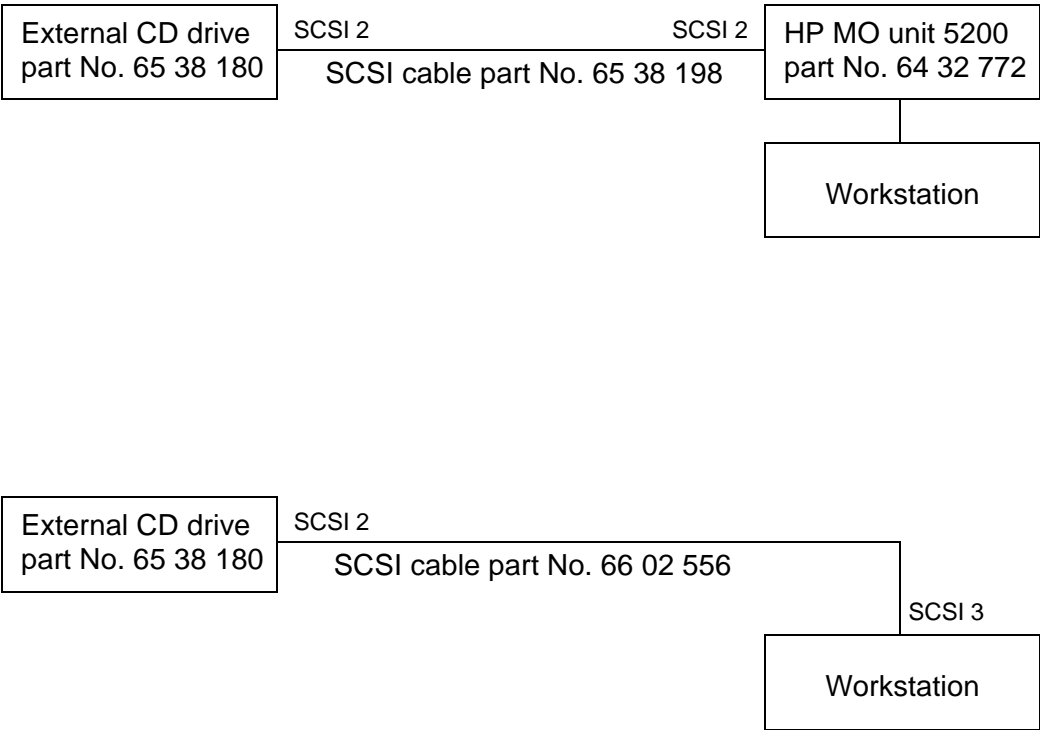


Fig. 2 SCSI cables

Database log file

This is what a log file (see Page 5 - 28) might look like.

```
=====
Checking patients
=====
-----
Reading patient id: 00-00001
Last name: Mrs A
First name: A
Date of birth: 400101
Gender: female
Number of examinations: 0          NOTICE! Not used (always 0).
-----
Reading patient id: 00-00002
Last name: Mr B
First name: B
Date of birth: 500101
Gender: male
Number of examinations: 0
-----
Reading patient id: 00-00003
Last name: Mrs C
First name: C
Date of birth: 600101
Gender: female
Number of examinations: 0
-----

=====
Checking folders
=====
-----
Reading folder id: 877524428
patient id: 00-00001
Folder name:
Folder type: spot examination
Folder creation time: 877524428
Folder creator: user
Examination type: Spot
Examination subject: Unknown
Number of images: 0          NOTICE! This folder is empty.
Stereo data: 0
Target side: SIN
-----
Reading folder id: 877526913
patient id: 00-00002
Folder name:
Folder type: spot examination
Folder creation time: 877526913
```

Folder creator: user
 Examination type: Spot
 Examination subject: Unknown
 Number of images: 1
 Stereo data: 0
 Target side: SIN

 Reading folder id: 877527322
 patient id: 00-00003
 Folder name:
 Folder type: spot examination
 Folder creation time: 877527322
 Folder creator: user
 Examination type: Spot
 Examination subject: Unknown
 Number of images: 3
 Stereo data: 0
 Target side: SIN

 Reading folder id: 877567455
 patient id: 00-00003
 Folder name:
 Folder type: stereo examination
 Folder creation time: 877567455
 Folder creator: user
 Examination type: Core biopsy
 Examination subject: Calcifications
 Number of images: 3
 Stereo data: 0
 Target side: SIN

 Reading folder id: 982139514
 patient id:
 Folder name: Stereo Calibration
 Folder type: stereo examination
 Folder creation time: 982139514
 Folder creator: user
 Examination type: Stereo Calibration
 Examination subject:
 Number of images: 4
 Stereo data: 0
 Target side: SIN

NOTICE! Empty patient id, used for biopsy calibration.

 =====
 Checking filestore
 =====

 Reading folder id: 877524428
 Copied to MO disk #0
 Last modified 877524428

NOTICE! There are no images to save,
 see Page A2-1.

 Reading folder id: 877526913

```
Copied to MO disk #0
Last modified 877526913
-----
Reading folder id: 877527322
Copied to MO disk #2
Last modified 877527322
-----
Reading folder id: 877567455
Copied to MO disk #3
Last modified 877567455
-----
Reading folder id: 982139514
Copied to MO disk #3
Last modified 982140155
-----
=====
Checking patientfolder
=====
-----
Reading patient id:
Number of folders: 1
Folder id 0:          982139514
-----
Reading patient id: 00-00001
Number of folders: 1
Folder id 0:          877524428
-----
Reading patient id: 00-00002
Number of folders: 1
Folder id 0:          877526913
-----
Reading patient id: 00-00003
Number of folders: 2
Folder id 0:          877527322
Folder id 1:          877567455
-----
```

NOTICE! Images not saved to local MO disk yet,
see Page A2-2.

NOTICE! Empty patient id, used for biopsy
calibration.

This page intentionally left blank.

Customer specific data

Parameter	Value
Language	
ID Pattern	
Date Pattern	
Service Center	
Fictive patient ID	
Institution Name	
Department Name	

Test protocol repetitiveness**Biopsy calculations**

Tester	Object	X	Y	Z	ΔX	ΔY	ΔZ
	1						
	1						
	1						
	1						
	2						
	2						
	2						
	2						
	3						
	3						
	3						
	3						
	4						
	4						
	4						
	4						

Test protocol CCD camera calibration

CSE in charge: _____

Calibration date: _____

TableGrid ☐Non grid ☐**For grid table only**

Grid fast speed time _____ (1500 ms)

Grid fast speed _____ (99% of max)

Grid slow speed _____ (40% of max)

For grid or non grid table

PMMA _____ (40 mm grid, 30 mm non grid)

kV _____ (27 kV grid, 26 kV non grid)

Anode/filter material _____ (Mo/Mo grid, Mo/Mo non grid)

Mode _____ (AEC mode grid, AEC mode non grid)

Test protocol image quality

AEC function

Grid table

Measured mean value, using 20 mm PMMA

Anode/Filter	26 kV	28 kV	30 kV	32 kV
Mo/Mo				

Measured mean value, using 40 mm PMMA

Anode/Filter	26 kV	28 kV	30 kV	32 kV
Mo/Mo				

Measured mean value, using 60 mm PMMA

Anode/Filter	26 kV	28 kV	30 kV	32 kV
Mo/Mo				
Mo/Rh				
W/Rh				

Measured mean value, using 80 mm PMMA

Anode/Filter	26 kV	28 kV	30 kV	32 kV
Mo/Rh				
W/Rh				

Non grid table

Measured mean value, using 20 mm PMMA

Anode/Filter	26 kV	28 kV	30 kV	32 kV
Mo/Mo				

Measured mean value, using 40 mm PMMA

Anode/Filter	26 kV	28 kV	30 kV	32 kV
Mo/Mo				

Measured mean value, using 60 mm PMMA

Anode/Filter	26 kV	28 kV	30 kV	32 kV
Mo/Mo				
Mo/Rh				
W/Rh				

Measured mean value, using 80 mm PMMA

Anode/Filter	26 kV	28 kV	30 kV	32 kV
Mo/Rh				
W/Rh				

Resolution**Normal resolution**

kVp used during tests

 kVp

mAs measured during test exposure with AEC and 45 mm PMMA

 mAs

mAs used during consecutive tests

 mAs

Anode/filter combination: Mo/Mo

Resolution: Normal

Best resolution visible lp/mm*Performance criteria:* ≥ 10 lp/mm

Anode/filter combination: W/Rh

Resolution: Normal

Best resolution visible lp/mm*Performance criteria:* ≥ 10 lp/mm

Comments: _____

High resolution (if used clinically)

kVp used during tests

 kVp

mAs measured during test exposure with AEC and 45 mm PMMA

 mAs

mAs used during consecutive tests

 mAs

Anode/filter combination: Mo/Mo

Resolution: High

Best resolution visible lp/mm*Performance criteria:* ≥ 13 lp/mm

Anode/filter combination: W/Rh

Resolution: High

Best resolution visible lp/mm*Performance criteria:* ≥ 13 lp/mm

Comments: _____

Test protocol protective earth measurement

Instrument type: _____

Serial No.: _____

Calibration date: _____

Biopsy unit

Measuring point	Measured resistance (Ω)	Resistance limit (Ω)
1		0.1
2		0.1
3		no limit
4		0.1
5		0.1
6		no limit
7		0.1

Biopsy controller

Measuring point	Max. measured resistance (Ω)	Resistance limit (Ω)
All		0.1

CCD camera

Measuring point	Max. measured resistance (Ω)	Resistance limit (Ω)
All		0.1

Biopsy controller cable duct

Measuring point	Max. measured resistance (Ω)	Resistance limit (Ω)
All		0.1